OPTIKA

M I C R O S C O P E S
I T A L Y

Ver. 3.0.1

B-380 - B-380ALC

OPERATION MANUAL

GUIDA UTENTE

MANUEL D'INSTRUCTIONS

MANUAL DE INSTRUCCIONES

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

OPTIKA MICROSCOPES - ITALY

www.optikamicroscopes.com - info@optikamicroscopes.com

1.0 DESCRIPTION	pag. 3
2.0 INTRODUCTION	pag. 5
3.0 UNPACKING AND ASSEMBLY	pag. 6
4.0 USING THE MICROSCOPE	pag. 7
4.1 Observation head	
4.2 Place the specimen on the stage	
4.3 Illumination system settings	
4.4 Interpupillary distance	
4.5 Focus tension adjustment	
4.6 Diopter adjustment	
4.7 Condenser	
4.8 Numerical aperture setting	
4.9 Phase rings centering (models B-380PH)	
5.0 AUTOMATIC LIGHT CONTROL (B-382-ALC MODELS)	pag. 10
6.0 MAINTENANCE	pag. 11
6.1 Microscopy environment	
6.2 Before and after using the microscope	
6.3 Precautions for a safe use	
6.4 Cleaning the optics	
7.0 RECOVERY AND RECYCLING	pag. 12

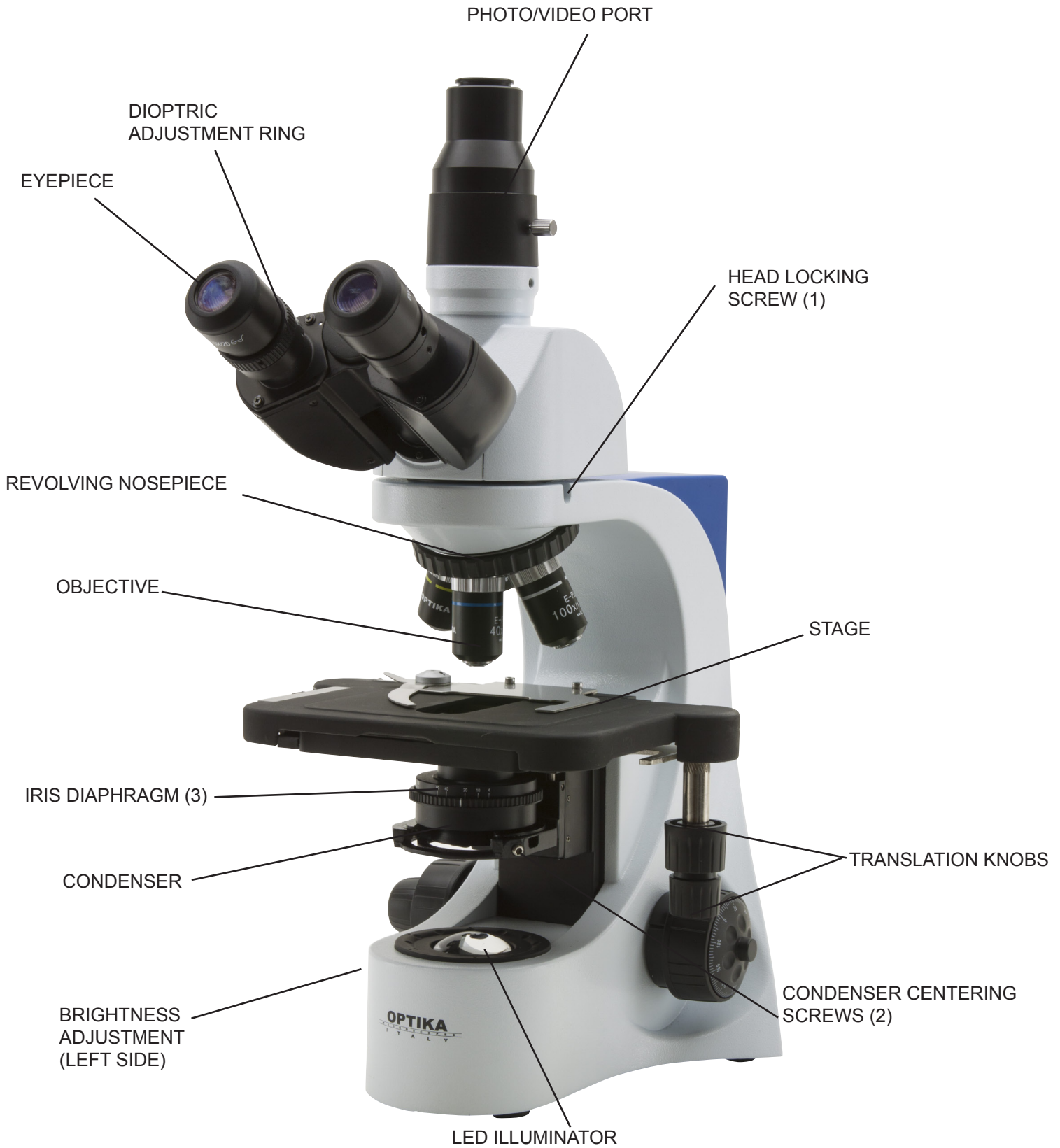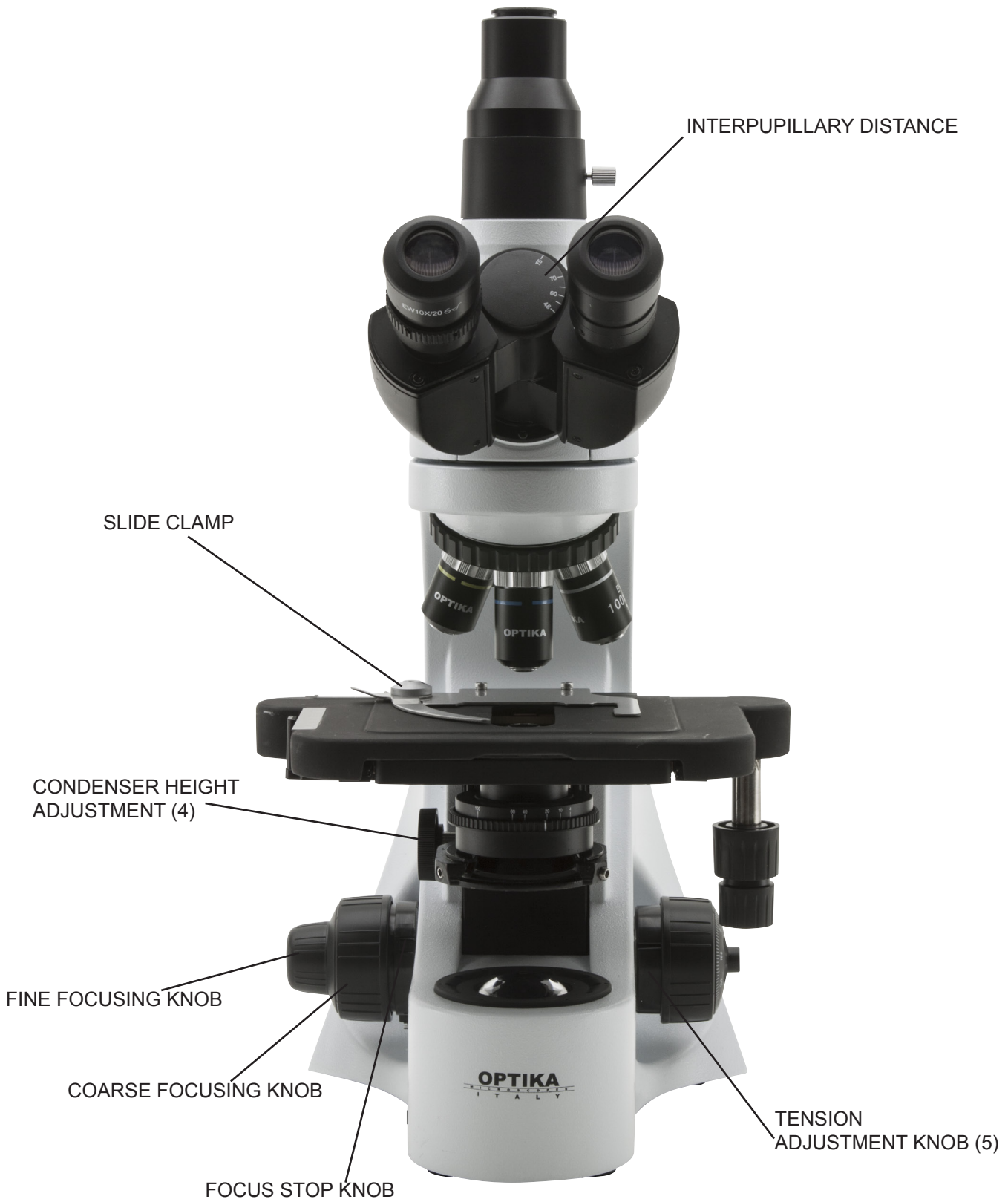This microscope is a scientific precision instrument designed to last for many years with a minimum of maintenance. It is built to high optical and mechanical standards and to withstand daily use.

Optika reminds you that this manual contains important information on safety and maintenance, and that it must therefore be made accessible to the instrument users.

Optika declines any responsibility deriving from instrument uses that do not comply with this manual.

2.1 Safety guidelines

This manual contains important information and warnings regarding safety about installation, use and maintenance of the microscope B-380. Please read this manual carefully before using the equipment. To ensure safe use, the user must read and follow all instructions in this manual. OPTIKA products are designed for safe use in normal operating conditions. The equipment and accessories described in the manual are manufactured and tested according to industry standards for safety instrumentation laboratory. Misuse can cause personal injury or damage to the instrument. Keep this manual at hand close to the instrument, for an easy consultation.

2.2 Electrical safety

Before connecting the power cord to wall outlet, ensure that your mains voltage for your region corresponds to the voltage supply of the instrument, and that the illuminator's switch is in position OFF. The user must observe the safety regulations in force in his region. The instrument is equipped with CE safety marking, in any case the user has full responsibility concerning the safe use of that instrument.

2.3 Warning/Caution symbols used in this manual

The user should be aware of safety aspects when using the instrument. Warning or hazard symbols are shown below. These symbols are used in this manual.

DANGER

The instructions on this symbol to avoid possible severe personal injuries.

WARNING

Warning of use; the incorrect operation on the instrument can cause damages to the person or instrument.

WARNING

Possibility of electric shock.

HOT!

Attention: high temperature surfaces. Avoid direct contact.

NOTE

Technical notes or usage tips.

3.0 UNPACKING AND ASSEMBLY

The microscope is located in a styrofoam moulded packaging. After removing the adhesive tape from all packaging, lift the top half of the packaging. Pay attention not to drop or damage the optical components (objectives and eyepieces). Extract the microscope from its packaging with both hands (one around the top arm and one around the base) and place it on a stable surface. Keep it away from solvents, chemical vapors and excessive moisture. Avoid high temperature environments, the direct sunshine and excessive vibrations, which could affect the performance of the instrument.

3.1 Operating environment

Temperature: : 10 - 36°C (50 – 96.8°F)

Relative humidity: 0 – 85% up to 30°C (86°F)

3.2 Unpacking microscope

Control the packaging to ensure that all material is present. We recommend that you take note of all the accessories to facilitate any future orders of spare parts and technical support calls. Make sure that in the packaging no small accessories or small parts remain. Please keep the original packaging in a safe place for future transport needs of microscope or accessories.


Never touch the glass surfaces such as lenses or filters. Traces of grease or other residues can reduce the vision quality of the final image and corrode the surface of lenses in a short time.

3.3 Installing the microscope

Set the optical head on the top arm through the locking screw. Insert the eyepieces into the tubes and lock them with the small screws which are located to the side of the tubes. Remove the protective film from the stage of the microscope.

3.4 Connect the external power supply into the socket at the base

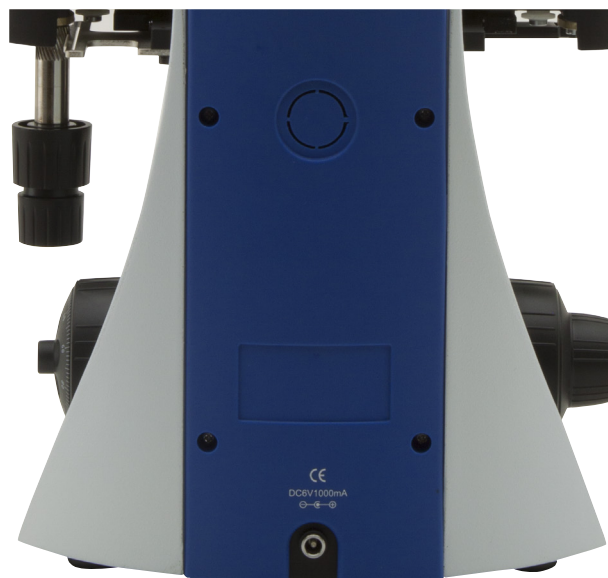Once positioned and installed with the necessary components, the microscope is ready to be used. Your microscope is a laboratory instrument designed to last. Handle it always carefully and avoid abrupt vibrations or shocks.

Always disconnect the power cable from the microscope when not in use for long time, while you clean it or when you perform any maintenance.

WARNING

AVOID DISASSEMBLING THE INSTRUMENT

Do not disassemble the instrument. This entails the cancellation of the warranty and may cause malfunction.

4.1 Observation head

Loosen the lock-screw (1), turn the observation head to a comfortable position for observation, and then lock the lock-screw.

4.2 Place the specimen on the stage

Lock the specimen slide on the mechanical stage using the slide clamp. Ensure that the specimen is centred over the stage opening by adjusting the coaxial knobs of the stage.

4.3 Illumination system settings

The microscope is fitted with a white LED illuminator. Before turning on the illumination system, read the section 5.3 about electrical safety precautions. Insert the plug of the cable into the power socket and turn on the switch on the side of the main body. Turn the brightness adjustment knob to a brightness suitable for observation.

4.4 Adjust interpupillary distance

Hold the right and left parts of the observation head with both hands and adjust the interpupillary distance by turning the two parts until one circle of light can be seen.

The white dot (°) placed on the left eyepiece shows the set interpupillary distance. Just remember this value to help on later settings.

4.5 Focus tension adjustment

The tension of the coarse focusing knob is preset by factory. To change the tension according to your preference, just rotate the knob (5) clockwise in order to increase it.

Excessive tension could damage the mechanism of focus. A too loosed tension causes the descent of the stage by gravity or a sudden loss of focus. In this case, rotate the knob (5) to increase the tension.

4.6 Diopter adjustment


Turn the dioptic adjustment ring on the left eyepiece to the zero position. Turn the coarse focus knob in order to focus the slide with an objective with low magnification. Adjust the fine focus knob until you obtain a clear and defined picture observing with the right eye, and then act on the left dioptic compensation ring observing with the left eye. When the image appears in focus, choose the necessary objective with the revolving nosepiece.

4.7 Condenser

Raise or lower the condenser through the knob (4) to obtain a clear and uniform illumination of the sample.

To center the condenser: completely close the iris diaphragm (3). Using the condenser centering screws (2), move the diaphragm in the center of the field of view. Then gradually expand the diaphragm until it is tangent to the edges of the field of view. If necessary, you can perform an additional adjustment.

The condenser is centered when the edges of iris diaphragm are tangent to the field of view.

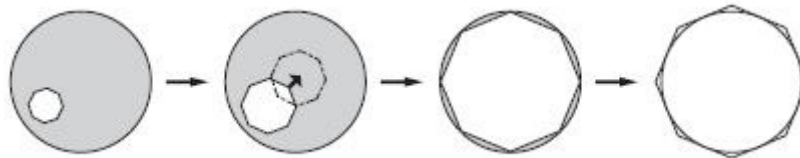4.8 Numerical aperture setting

The value of the numerical aperture (N.A.) of the diaphragm is an indication of the contrast of the illumination system. Matching the value of illumination system's N.A. with that of the objective ensures the best results in terms of contrast and image quality. To set the numerical aperture of the Illuminator, adjust the opening of the iris diaphragm (3). In this way you control contrast and image resolution. For samples with low contrast set the iris to about 75% of the value of the objective's numerical aperture.

4.9 Phase rings centering (models B-380PH)

For models equipped with phase contrast set, you have to perform the centering of the phase rings. Remove an eyepiece from the head and insert the centering telescope in the empty tube.

Insert the 10x objective rotating the nosepiece.


Rotate the turret of the condenser until you reach the inscription "10".

Loosening the lock screw of the centering telescope, focus on the light ring that you observe.

Rotate the two centering screws in order to move the bright ring until it is perfectly aligned with the dark ring.

Repeat for the other objectives (only as a verification of the correct centering).

Phase rings will be centered when you see an image like this:

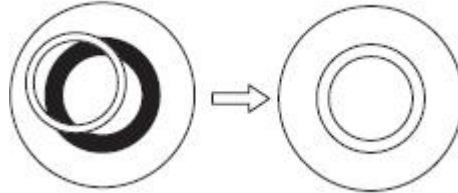B-382-ALC models have a special function called “Automatic Light Control”.

The level of light is adjusted by the microscope in order to maintain the same level as the one the user has chosen, no matter if the aperture of the diaphragm changes, another objective is inserted, opacity of the sample changes, etc.

1) Set the focus on the sample as described in the previous chapters, using the objective of your choice.

2) Rotate the brightness adjustment knob in order to get a comfortable level of illumination. NOTE: ALC system works properly with medium-to-high light intensity, do not set the brightness to a minimum level.

3) Press the ALC button on the side of the microscope.

Now the ALC system is working: if you open or close the aperture diaphragm or change objective, the system will try to maintain the same level of illumination as the one you stored when you pressed the button.

4) Press the ALC button again to return to a manual control of illumination.

NOTE: ALC system is not suitable for use with phase contrast.

6.1 Microscopy environment

This microscope is recommended to be used in a clean, dry and shock free environment with a temperature of 0-40°C and a maximum relative humidity of 85 % (non condensing). Use a dehumidifier if needed.

6.2 Before and after using the microscope

- The microscope should always be kept vertically when moving it and be careful so that no moving parts, such as the eyepieces, fall out.
- Never mishandle or impose unnecessary force on the microscope.
- Never attempt to service the microscope yourself.
- After use, turn off the light immediately, cover the microscope with the included dust-cover, and keep it in a dry and clean place.

6.3 Precautions for a safe use

- Before plugging in the power cord with the supply, make sure that the supplying voltage of your region matches with the operation voltage of the equipment and that the lamp switch is in off-position.
- Users should observe all safety regulations of the region. The equipment has acquired the CE safety label. However, users do have full responsibility to use this equipment safely.

6.4 Cleaning the optics

- If the optical parts need to be cleaned try first to: use compressed air.
- If that is not sufficient: use a soft lint-free piece of cloth with water and a mild detergent.
- And as a final option: use the piece of cloth moistened with a 3:7 mixture of ethanol and ether.
Note: ethanol and ether are highly flammable liquids. Do not use them near a heat source, near sparks or near electric equipment. Use these chemicals in a well ventilated room.
- Remember to never wipe the surface of any optical items with your hands. Fingerprints can damage the optics.
- Do not disassemble objectives or eyepieces in attempt to clean them.

If you need to send the microscope to Optika for maintenance, please use the original packaging.

Art.13 Dlsg 25 July 2005 N°151. “According to directives 2002/95/EC, 2002/96/EC and 2003/108/EC relating to the reduction in the use of hazardous substances in electrical and electronic equipment and waste disposal.”

The basket symbol on equipment or on its box indicates that the product at the end of its useful life should be collected separately from other waste.

The separate collection of this equipment at the end of its lifetime is organized and managed by the producer. The user will have to contact the manufacturer and follow the rules that he adopted for end-of-life equipment collection. The collection of the equipment for recycling, treatment and environmentally compatible disposal, helps to prevent possible adverse effects on the environment and health and promotes reuse and/or recycling of materials of the equipment. Improper disposal of the product involves the application of administrative penalties as provided by the laws in force.

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