

## Jaw Crusher BB 200

### General Information

The Jaw Crusher BB 200 is used for the rapid, gentle crushing and pre-crushing of medium-hard, hard, brittle and tough materials. Its variety of materials offered, including heavy-metal free steel, its efficiency and safety makes it ideal for sample preparation in laboratories and industrial plants.

For small amounts of sample the BB 200 can be used batch-wise; for larger amounts it can be operated continuously. Control of the gap width and zero point adjustment allow for reproducible results. The Belleville spring washer protects the unit from overloads and in conjunction with the central lubrication guarantees a long service life.

### Application Examples

alloys, basalt, cement clinker, ceramics, chamotte, coal, coke, construction materials, feldspar, glass, granite, minerals, ores, oxide ceramics, quartz, rocks, silicon, slag, ...

### Product Advantages

- excellent crushing performance
- connector for dust extraction
- central lubrication
- suitable for integration in automatic installations
- special version for grinding semi-conductor materials
- wear compensation with zero-point adjustment
- gap width setting
- overload protection
- wide range of materials for contamination free grinding
- no-rebound feed hopper with quick-release clamp
- brake motor with safety switch
- easy-to-clean crushing chamber
- continuous grinding



## Jaw Crusher BB 200

### Features


Applications	coarse and pre-crushing
Field of application	chemistry / plastics, construction materials, engineering / electronics, environment, geology / metallurgy, glass / ceramics
Feed material	medium-hard, hard, brittle, tough
Size reduction principle	pressure
Material feed size*	< 90 mm
Final fineness*	< 2 mm
Material of grinding tools	manganese steel, stainless steel, tungsten carbide, steel 1.1750 (for heavy-metal free grinding)
Jaw width	90 x 90 mm
Gap width setting	0 - 30 mm
Gap width display	scale
Zero point adjustment	yes
Hinged hopper	yes
Dust extraction unit	yes
Central lubrication	optional
Process line version	optional
Collector capacity	5 l
Drive	1-phase motor / 3-phase motor
Drive power	1.5 kW
Protection code	IP 54
W x H x D closed	450 x 1160 x 900 mm
Net weight	~ 300 kg
Workplace related emission value	LpAeq 92 dB(A)
Documentation	Operation & Application Video
Standards	CE

\*depending on feed material and instrument configuration/settings

### Function Principle

The Jaw Crusher BB 200 is a robust and powerful forced-feed crusher. The feed material passes through the no-rebound hopper and enters the crushing chamber. Size reduction takes place in the wedged-shaped area between the fixed crushing arm and one moved by an eccentric drive shaft. The elliptical motion crushes the sample which then falls under gravity.

As soon as the sample is smaller than the discharge gap width, it falls into a removable collector. The continuous gap width setting with scale ensure optimal size reduction in accordance with the set gap width value

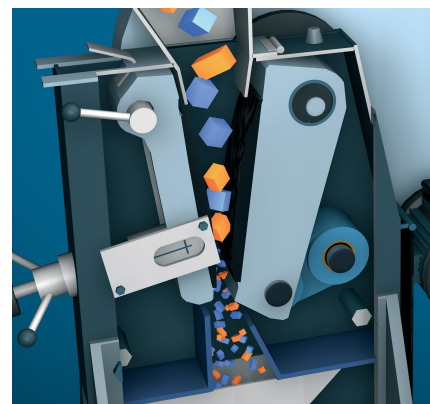



**Wolflabs**

# Wolf Laboratories Limited

[www.wolflabs.co.uk](http://www.wolflabs.co.uk)

Tel: 01759 301142

Fax: 01759 301143

[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)



**Use the above details to contact us if this literature doesn't answer all your questions.**

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

