



EL - LT / EL - LTS Low Temperature



The widest range in low temperature freezers

# The widest range in low temperature freezers



The LT comes with a handle lock



The pull-out defrost drain makes defrosting easy



# Temperature freezers



In laboratories and for storing special types of ice cream, low temperatures are often essential. And if you want to preserve fish for e.g. sushi, stable low temperatures are the best guarantee for keeping a high quality.

The LT/LTS-range from Elcold offers temperatures down to  $-45^{\circ}$  and freezers from 100 to almost 500 litres volume.

The temperature is kept at a constant level and can easily be monitored on the digital display.

Setting and programming is also done on the controller, which furthermore offers the possibility for adding an external alarm.

The insulation of 100 mm polyurethane foam ensures that the temperature will only rise very slowly, in case of power failure.

Together with the traditional high Elcold quality, that is the best insurance you can get.

The electronic thermostat is easy to read and makes it easy to monitor the freezer's status




The LTS is delivered with glass lids



## Technical data

		<b>EL 11 LT/LTS</b>	<b>EL 21 LT/LTS</b>	<b>EL 31 LT/LTS</b>	<b>EL 41 LT/LTS</b>	<b>EL 51 LT/LTS</b>
Gross volume	litres	130/140	226/243	299/322	358/384	416/447
Net volume	litres	130/116	226/203	299/270	358/323	416/356
Gross volume	cu.ft.	4,6/5,0	8,1/8,7	10,7/11,5	12,8/13,7	14,9/16,0
Net volume	cu.ft.	4,6/4,1	8,1/7,3	10,7/9,6	12,8/11,5	14,9/12,7
<b>A</b>	mm	725	1055	1305	1505	1705
<b>B</b>	mm	655	655	655	655	655
<b>C</b>	mm	865	865	865	865	865
<b>D</b>	mm	780	780	780	780	780
<b>E</b>	mm	520	850	1100	1300	1500
<b>F</b>	mm	320	650	900	1100	1300
<b>G</b>	mm	450	450	450	450	450
<b>H</b>	mm	650/695	650/695	650/695	650/695	650/695
<b>J</b>	mm	400/445	400/445	400/445	400/445	400/445
<b>L</b>	mm	850	850	850	850	850
Isulation thickness	mm	100	100	100	100	100
Climate class 4	(EN 441)	+	+	+	+	+
Handle with lock	(LT)	+	+	+	+	+
Defrost drain		+	+	+	+	+
Electronic thermostat		+	+	+	+	+
Acoustic alarm		+	+	+	+	+
Temp. range at 30°C amb. °C	(LT)	-10° - -45°	-10° - -45°	-10° - -45°	-10° - -45°	-10° - -45°
Temp. range at 30°C amb. °C	(LTS)	-10° - -30°	-10° - -30°	-10° - -30°	-10° - -30°	-10° - -30°

Subject to alterations. The technical specifications should be taken as a guidance.

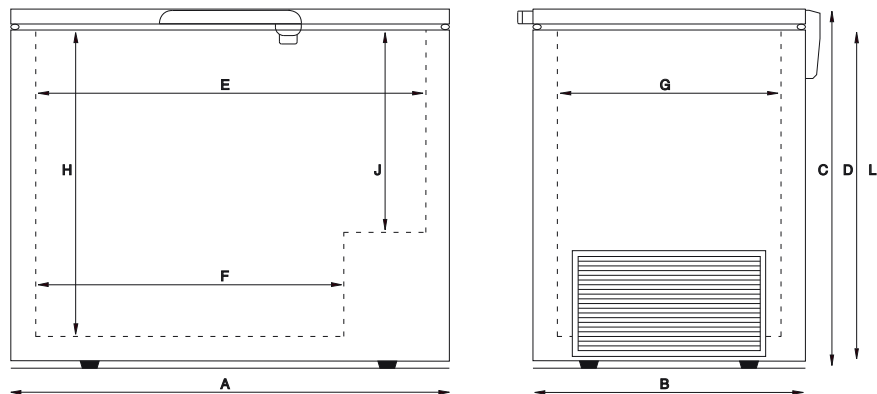

Wolf Laboratories Limited, Colenso House, 1 Deans Lane, Pocklington, York, YO42 2PX

Phone: 01759 301142 Fax: 01759 301143

Email: sales@wolflabs.co.uk www.wolflabs.co.uk