

# GX Wolf Microscopes

## XTX-7CW Series

## product datasheet



High Quality  
Stereo  
Microscope  
&  
Gemological  
Accessories

### Features and Specifications

#### Ergonomic

Binocular head which is ergonomically angled at 45 degrees for maximum comfort during prolonged observations.

The binocular head can be rotated 360 degrees.

To eliminate eyestrain each user can set up perfect viewing by adjusting the diopter setting on the left eyepiece tube.

The two ocular tubes can be adjusted to accommodate the perfect interpupillary distance for each user.

#### Great Design

The microscope head is mounted on sturdy and stylish microscope 'arm-type' metal stand.

#### Built to Last

The focus knob moves the head smoothly using a metal rack and pinion mechanism constructed from tough metal bearing for a long working life.

#### Superb Optics

The optical parts of the microscope consist of:

a pair of high quality 10X WF eyepieces. The XTX7CWKit version includes sets of additional eyepieces providing a very wide range of magnifications.

#### 3D Viewing

A rotatable objective nosepiece which contains 2 pairs of objectives. Each objective in each pair links directly to an eyepiece. The objectives are angled. This gives the user a great stereo 3D image of the specimen under observation. The optics are all made from high quality glass and produce crisp, bright images with minimal chromatic aberration.

#### 80X Magnification

The objective pairs have magnification factors of 2X and 4X to produce a total magnification of 20X and 40X when used with the 10X eyepieces and up to 80X when used with the 20X eyepieces.

#### Dual Illuminated Specimens

Illumination is provided by lamps above and below the specimen.

#### Independent Intensity Control

The lamps can be turned on independently and have variable intensity controls.

The base light is covered by a frosted glass plate providing a evenly illuminated specimen.


#### Two Versions with Full Gemological Kit

The XTX7CWKit version includes a set of gemological tools:

- o Jewel Tweezers which are mounted on a rotatable pillar which sits in the slide clip holes
- o Darkfield stage which is just perfect for observing gems and diamonds giving a black background while illuminating the gem from the sides. This stage sits on top of the base illuminator.

#### Spares and Warranty

Every microscope is supplied with spare bulbs and comes with a 1 year warranty.

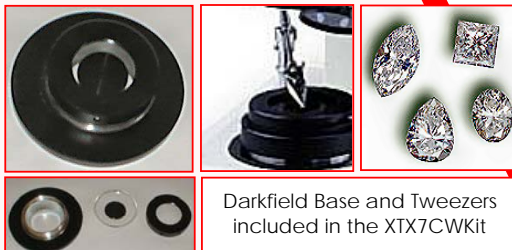

Darkfield Base and Tweezers included in the XTX7CWKit

Model:	XTX-7C-W	XTX-7C-WKit
Head type	Binocular	Binocular
Eyepieces	WF 10X	WF 10X (also see accessories)
Diopter Adjustment	Std	Std
Stand type	Arm	Arm
Illumination	Top & Base Adjustable	Top & Base Adjustable
Magnification with 10X eyepieces	20x & 40X	20x & 40X
<b>Accessories</b>		
Eyepieces (pair)	WF15x WF 20x Optional	WF15x WF 20x Standard
Gem Kit (darkfield base and jewellery tweezers)	Optional	Standard
Digital cameras	Optional	Optional

### Wolf Laboratories Limited

Colenso House  
1 Deans Lane  
Pocklington  
York YO42 2PX

Tel 01759 301142  
Fax 01759 301143

Email [sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)  
Web [www.wolflabs.co.uk](http://www.wolflabs.co.uk)

Due to our policy of continuous development specifications may change without notice

Jan 06