

THE **MINIMED90** X-RAY FILM PROCESSOR




FEATURES

- ◇ Compact size, ideal for small darkrooms
- ◇ Film sizes 10x10cm to 35x43cm
- ◇ 90 second processing cycle
- ◇ Energy and water conserving standby mode
- ◇ Fully automatic chemical replenishment and wash water control

Compact, versatile, cost effective film processor

SPECIFICATIONS

- ◇ Processing Time: 90 seconds (leading edge to leading edge)
120 seconds (dry to dry)
- ◇ Film Sizes: 10x10cm to 35x43cm
- ◇ Capacity: 90 films per hour (35x43cm)
- ◇ Film Types: Any film suitable for 90 second processing
- ◇ Chemistry: Any chemistry suitable for rapid processing
- ◇ Water consumption: 0.9 litres per minute (none during standby mode)
- ◇ Dimensions: Length: 83.8cm Width: 55.9cm Height: 62.2cm
- ◇ Weight: 45.5 kg

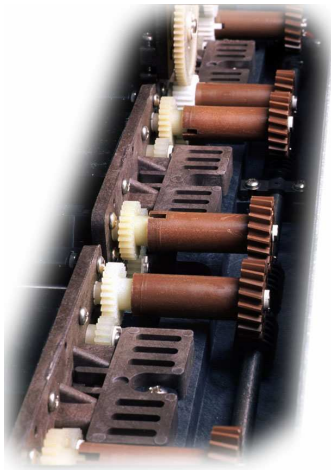

- ◇ Dual film sensors ensure all film sizes are detected
- ◇ Simple drive system for easy maintenance. No chain to stretch or break
- ◇ Front exit drier (standard) is easily removed for access to wash tank and for maintenance

OPTIONS

- ◇ Processor Stand (Height 74cm)
- ◇ Rear exit version for through-the-wall installation
- ◇ LED display panel interchangeable with on/off switch



LED panel is interchangeable with the on/off switch and indicates power and chemistry status

because your image matters