

## STAINING RACK ENDS AND RODS

Chrome plated brass ends with leveling feet for individual staining of 76mm slides over a sink. Stainless steel rods supplied in 1/4" diameter to lengths as required. Please specify in centimeters when ordering.


<b>Cat Ref</b>	
E96.6	Staining Rack Ends (Pair)
E96.62	Rods (cms)



## SLIDE WASHING TRAY

Stainless steel tray for washing sections in running water. Slides are placed at an angle and water passes through a number of holes in the base, flooding over the lip of the tray.

<b>Cat Ref</b>	<b>Size</b>
E89.06	263 x 300 x 62mm

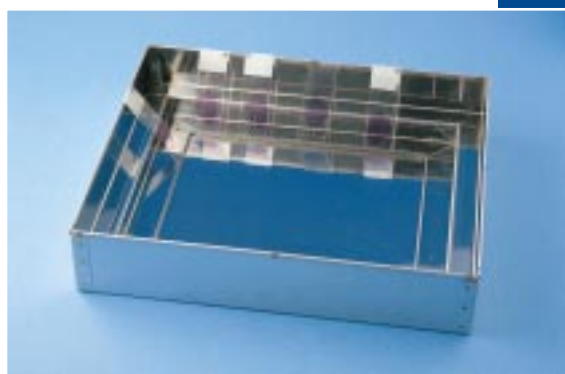


## IMMUNO SLIDE STAINING TRAYS

Developed for routine reagent techniques requiring a humid chamber, the apparatus allows simple, safe manipulation of 12 or 24 slides (76 x 26mm) using only one hand. Trays are made of stainless steel with a removable section supporting the slides. This section may be tilted to allow buffer or reagent to run off. The inner cover is of clear Perspex™ acrylic plate allowing visual examination and an extra cover of stainless steel prevents evaporation of reagents.

Cat Ref	Capacity	Dimensions (with lid)
E103.2	12 slides	350 x 140 x 35mm
E103.3	24 slides	350 x 255 x 35mm



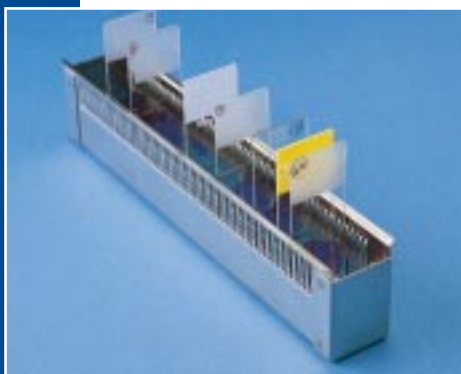

### VERTICAL SLIDE RACK

Stainless steel rack with grooves to hold 76 x 26mm microscope slides vertically.

Cat Ref	Dimensions	Capacity
E89	191 x 35 x 43mm	60 slides



### HORIZONTAL SLIDE RACK

Stainless steel rack with grooves to hold various microscope slides horizontally.


Cat Ref	Dimensions	Capacity
E89.05	237 x 80 x 24mm	60 slides



### SLIDE STAND

Stainless steel rack with four rails to support various microscope slides almost vertically.

Cat Ref	Dimensions	Capacity
E89.03	455 x 100mm	68 slides

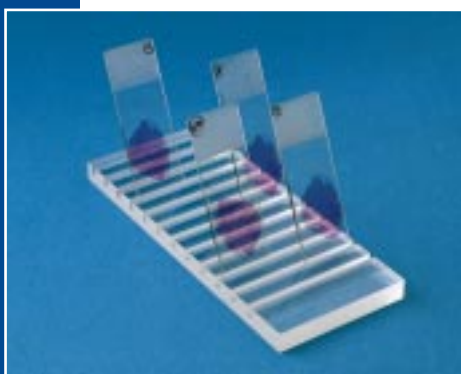

### PERSPEX™ ACRYLIC SLIDE RACK

Perspex™ acrylic plate with twelve transverse grooves to take microscope slides of various sizes.

Cat Ref	Dimensions	Capacity
E89.055	162 x 64mm	24 slides

NB – Capacities given are based on 76 x 26mm 1.0-1.2mm slides