


Helpful gadgets

For your Laboratory



Weighing dishes made of aluminium

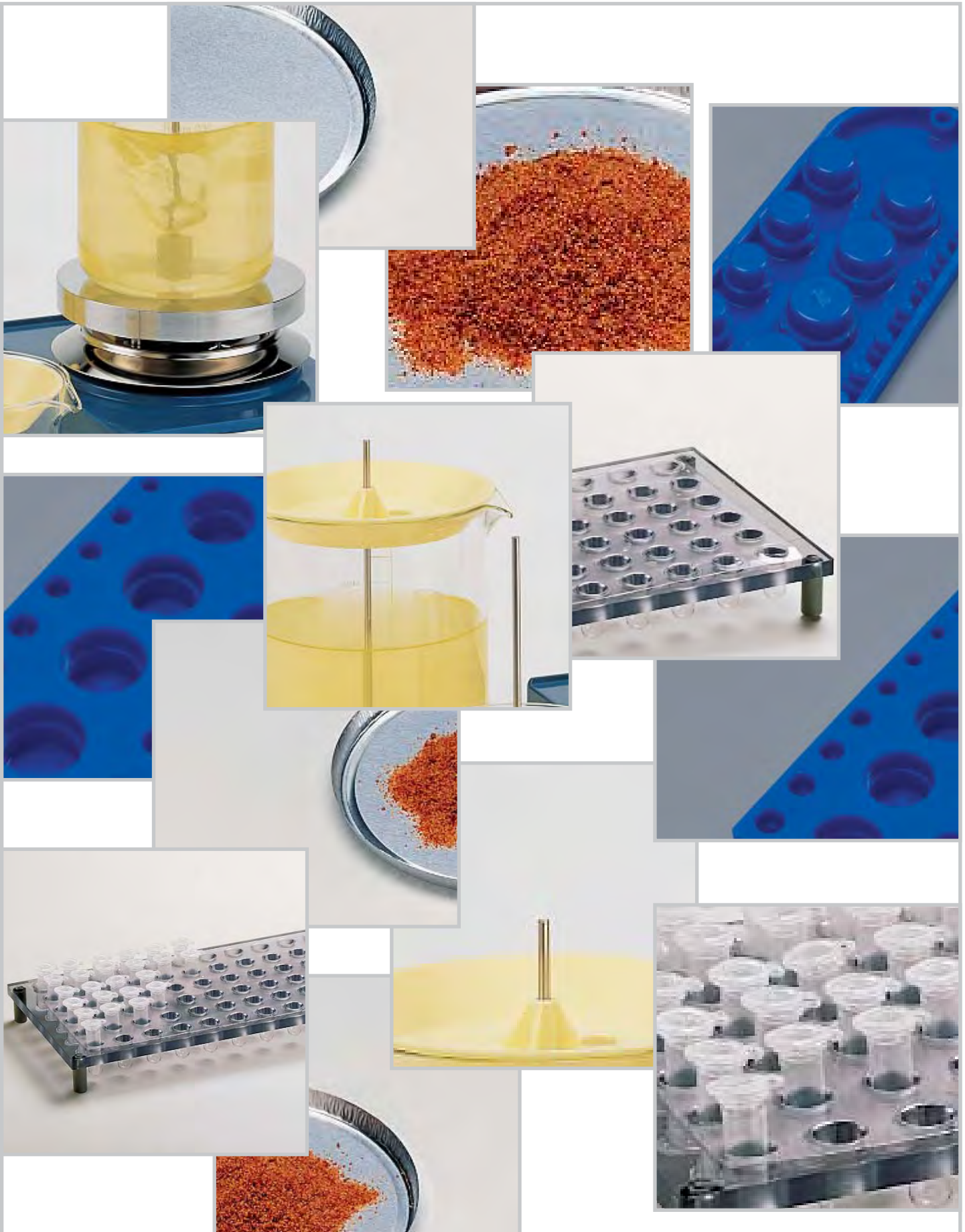

- weighing dishes, ideal for use in infrared dryers
- multifunctional

Ideal for moisture determination with infrared dryers. The samples should be distributed

horizontally and reflect the heat evenly over the complete dish area.

May also be used as laboratory dishes.

MAGNETO-STABIL Stirring inserts for beakers



- Smooth running, no spattering or evaporation of liquid when operating at high speed or continuously
- Minimal wear of the stirring bar, the beaker and the sample due to the device's high performance
- The conical lid safely prevents from spattering and reduces evaporation

Stirring inserts available in 3 different sizes suitable for all types of beakers from 100 ml up to 2000 ml capacity.

Designed with a chemical-resistant conical plastic cover with bushing and 2 pipetting bores as well as a teflon-coated magnetic stirring bar with mounting.

MESSFIX Tubing gauge



- small and handy
- easy measurement

For easy measurement of the outer or inner diameter of hoses, glass tubes etc. in the range of 1.5 to 16 mm diameter.


always ready at hand to check the diameter of a hose or tube.

To determine the inner diameter, the hose is simply shoved to a fitting spike and the diameter can be read.

A small and handy gauge made of polypropylene, that can easily be carried in every pocket. It is

To measure the outer diameter, the hose is stuck into a fitting hole and the diameter, too, can be read.

TESTIPLATT Test tube tray

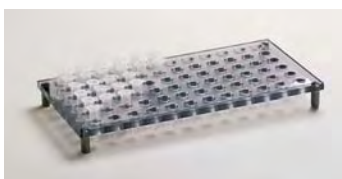

- ideal for serial assays
- resistant to acids and alkaline solutions

Made of hard PVC-foil, vacuum-formed. Ideal for serial assays (e. g. blood tests), resistant to acids and alkaline solutions. Transparent.

with round bottom, 10 mm diameter, 25 mm deep, 2 ml tube volume.

TESTIPLATT is very suitable as a support for reaction vessels.

Supporting trays made of transparent PVC add stability to TESTIPLATT.

72 (6x12) slightly conical tubes

Ordering Information

MAGNETO-STABIL 1
MAGNETO-STABIL 2
MAGNETO-STABIL 3

MESSFIX

TESTIPLATT

TESTIPLATT Supporting tray

Weighing dishes

for beakers 60 - 95 mm Ø, 200 mm height
for beakers 86 - 118 mm Ø, 200 mm height
for beakers 110 - 143 mm Ø, 200 mm height

120 x 12 x 50 mm (w x h x d)

275 x 25 x 137 mm (w x h x d)

100 mm Ø, 7 mm height (qty=80)

Cat.-No.

3.440 102

3.440 202

3.440 302

3.660 103

3.683 073

3.683 173

3.628 103

schuett[®]biotec.de

Rudolf-Wissell-Straße 13, D-37079 Göttingen, Germany
Fon +49 (0) 551/5 04 10-0, Fax +49 (0) 551/5 04 10-99
info@schuett-biotec.de
www.schuett-biotec.de