



*schuett* **petriturn**

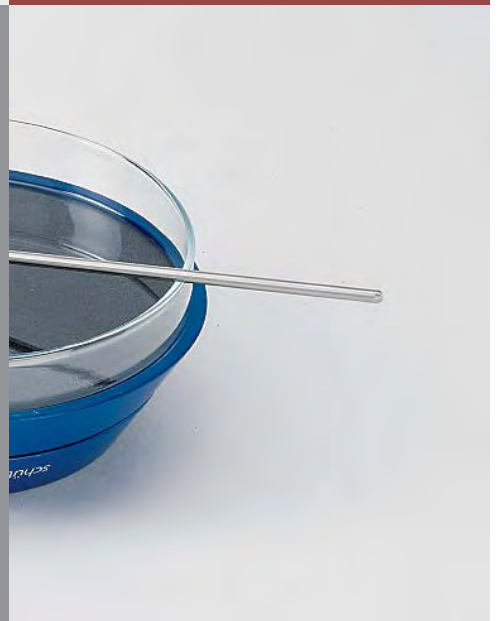

*Turntables for Petri dishes & Inoculating hooks*

*schuett* **biotec**.de

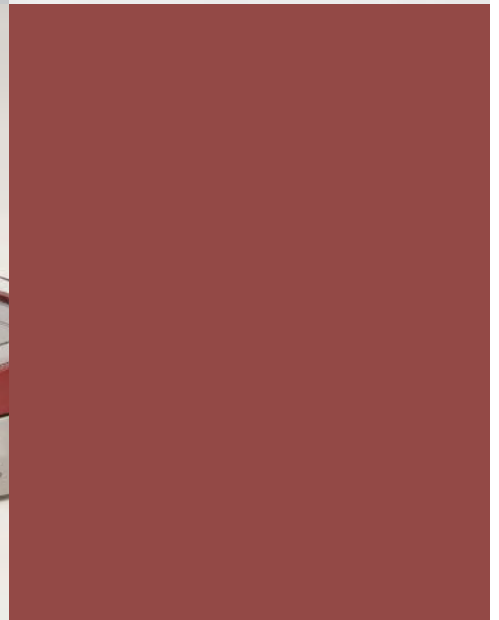

schuett petriturn-E with foot-switch and turntable-attachment for Petri dishes up to Ø 150 mm (optional)



Inoculating hook made of stainless steel (Drigalski-hook)



Inoculating hook made of glass (Drigalski-hook)

Easy & Uniform

### Petri dishes turntables

For uniform plating of liquids on agar-surfaces in Petri dishes. Easy to handle, tireless.

### schuett petriturn-E Electrically driven Petri dish turntable

With electrically drive and sensor. Body made of stainless-steel. Handling of the turntable contactless by sensor. For Petri dishes up to Ø 100 mm (up to Ø 150 mm optional).

Continuously adjustable, constant speed of 10-100 rpm and continuously adjustable operating time of 3 s to 2 min or permanent running.

Optional connection of a foot-switch.

Turntables for Petri dishes of other diameters available on request.

Comfortable

### schuett petriturn-M Manual operation

Powerful and long-running due to high mass of turntable.

Useable on both sides, that means two work spaces in one device by easiest turn-over. For Petri dishes with Ø 100 and Ø 150 mm diameter.



schuett petriturn-M with inoculating hook (optional)

### Technical data

#### petriturn-E

Dimensions (Ø x h)	160 x 70 mm
Weight	approx. 1 kg
Wall power supply	12 V DC 5 VA
Speed continuously, adjustable	up to 10-100 rpm
Operating time, continuously adjustable	3 s - 2 min / permanent running

#### petriturn-M

Dimensions (Ø x h)	160 x 45 mm
Weight	approx. 2.3 kg

### Ordering Information

schuett petriturn-E, electrically driven, 230 V	<b>Cat.-No.</b> 3.361 252
schuett petriturn-E, electrically driven, 115 V	3.361 262
schuett petriturn-M, manual operation	3.361 102

### Accessories

Foot-switch for schuett petriturn-E	3.361 302
Turntable-attachment for Petri dishes up to Ø 150 mm	3.361 402

### Inoculating hook

Inoculating hook made of stainless steel (Drigalski-hook), 180 mm length	3.362 152
Inoculating hook made of glass (Drigalski-hook), 145 mm length	3.362 162

**schuett**biotec.de

Rudolf-Wissell-Straße 13, D-37079 Göttingen, Germany  
Fon +49 (0) 551/5 04 10-0, Fax +49 (0) 551/5 04 10-99  
info@schuett-biotec.de  
www.schuett-biotec.de