


PORTABLE X-RAY UNIT



The optional Column and the Articulated Arm make this Portable X-ray Unit the ideal choice for portable / mobile radiography.

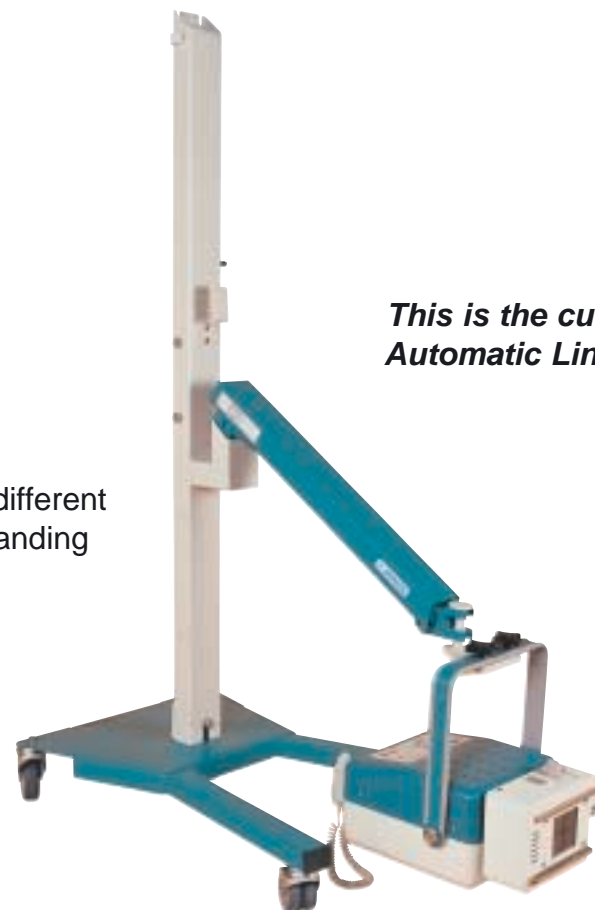

The Articulated Arm can be fixed at two different positions on the column, which allows standing (weight-bearing) exposures.



... up to 4 kW Power
110 kV

100 mA

300 kHz

0.6 mm

This is the currently the only unit on the market available with Automatic Line Compensation from 90 to 264 Volts (AC).



Comfortable and safe transportation in any working environment is ensured using the wheeled case with telescopic handle, available as an option with equipment.

The function of this high performance equipment is based on the following principle: for any combination of kV and mAs it always uses the **maximum mA value** with the **minimum Exposure Time**. In this way, the exposure parameters are automatically optimised.

- Exposure time is continuously System-calculated and displayed

The Unit automatically chooses the Focal Spot to apply based on the exposure parameters selected by the operator. Manual selection of Small and Large Focal Spot will be applied for those instances where the exposure values allow this.