



Wolflabs

Wolf Laboratories Limited

www.wolflabs.co.uk

Tel: 01759 301142

Fax:01759 301143

sales@wolflabs.co.uk




If the literature below doesn't answer all your questions, click these links to:

[Request details on an alternative product of the same type](#)

[Request details on accessories](#)

[Request details on other product categories](#)



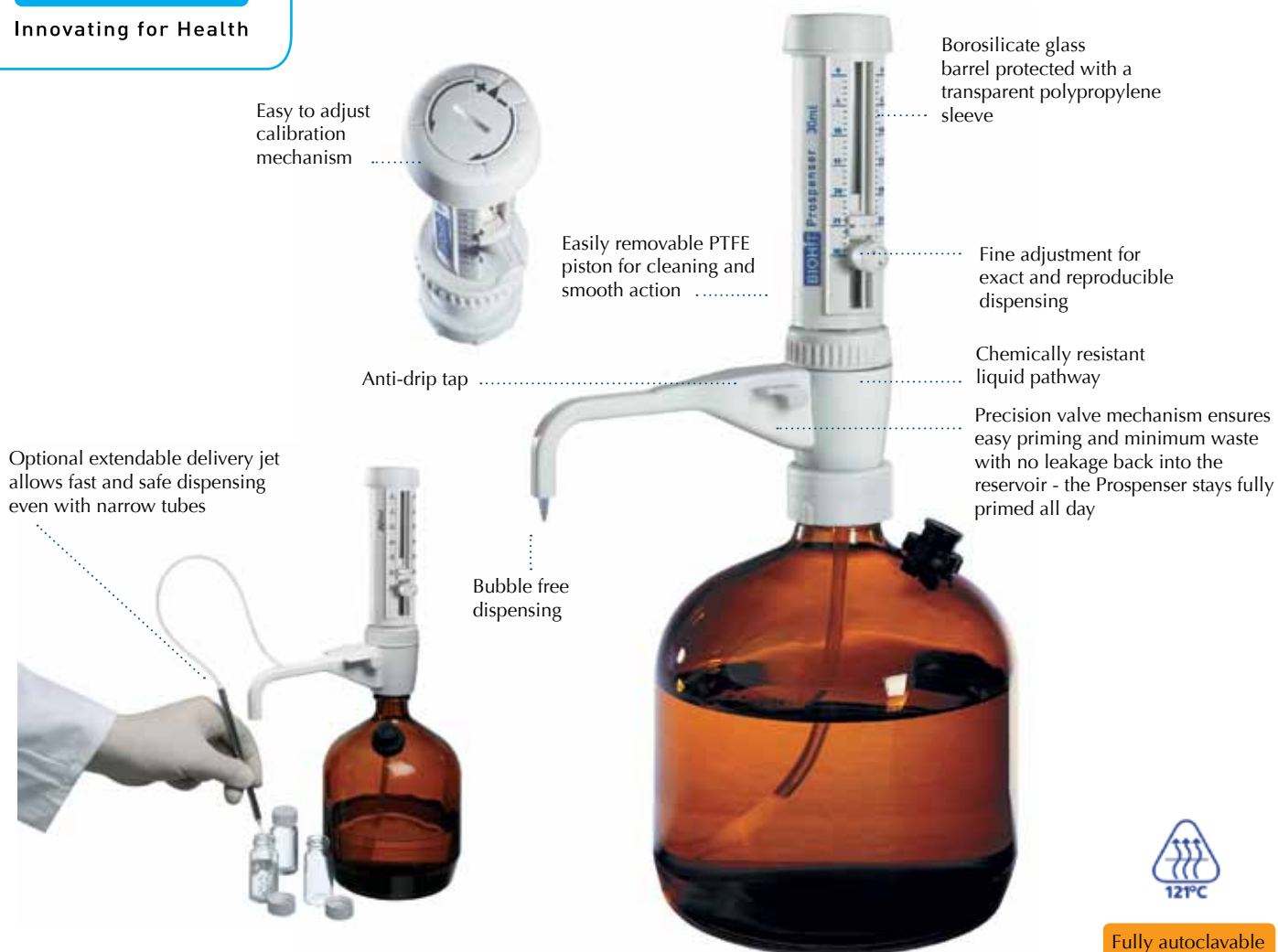

Prospenser Bottle-Top Dispenser

The easy-to-use Prospenser delivers trouble-free and reliable dispensing of liquids including strong acids, alkalis and solvents.

- All models are fitted with an anti-drip safety valve, a feature which guards against drips when the Prospenser is not in use.
- A fine adjustment mechanism allows for even greater reproducibility for repeat dispensing.
- Unlike other bottle-top dispensers, the Prospenser's glass barrel can be disassembled from the pedestal for thorough cleaning. This also ensures that in the unlikely event of

damage, you can replace the glass barrel as a spare part rather than have to buy a new dispenser.

- The head can be conveniently rotated 360° for operator safety and comfort.
- Wide range of adapters are included to fit the most common bottle sizes.
- Optional extendable delivery jet allows fast and safe dispensing even to narrow tubes.

Prospenser Performance specifications and ordering details

Cat. No.	Item	Increment	Max Volume	Inacc. (%)	Impr. (%)
723049	Prospenser 0.01–2.5 ml (with 38, 40 and 45 mm adaptors)	0.05 ml	2.5 ml	0.3 %	0.1 %
723050	Prospenser 0.1–5 ml (with 38, 40 and 45 mm adaptors)	0.1 ml	5 ml	0.3 %	0.1 %
723051	Prospenser 0.2–10 ml (with 38, 40 and 45 mm adaptors)	0.2 ml	10 ml	0.3 %	0.1 %
723052	Prospenser 1–30 ml (with 38, 40 and 45 mm adaptors)	1.0 ml	30 ml	0.3 %	0.1 %
723053	Prospenser 1–50 ml (with 38, 40 and 45mm adaptors)	1.0 ml	50 ml	0.3 %	0.1 %
721998	Extendable delivery jet (see picture)				

Please note that the Prospenser is not supplied with a reservoir





Biotrate Digital Burettes 0–30 and 0–50 ml

The streamlined Biotrate digital burette and dispenser delivers accurate, precise and convenient bottle-top titration, as well as optimum operator safety.

- The head can be conveniently rotated 360° for operator safety and comfort.
- See-through barrel ensures bubble-free titration.
- Automatic switch-off without losing the previous titrated volume, when left unattended
- No dripping when in static mode
- A new safety pedestal valve allows the Biotrate fluid path to be turned off manually to prevent accidental delivery. Its right-angled spout is also designed to enhance safety in use.
- Robust outer casing does not move during the priming or dispensing sequences.
- PTFE plunger and borosilicate glass barrel allow dispensing of aqueous solutions and some aggressive liquids.
- Powered by a 3.6V lithium battery giving a minimum life of 60,000 three-minute titrations. Automatic low battery voltage is indicated on the large easy-to-read LCD window.

Biotrate Performance specifications and ordering details

Cat. No.	Item	Increment	Inacc. (%)	Impr. (%)
723054	Biotrate 0–30 ml (with 33, 38 and 45mm adaptors)	0.01 ml	0.2 %	0.1 %
723055	Biotrate 0–50 ml (with 33, 38 and 45mm adaptors)	0.01 ml	0.2 %	0.1 %
721998	Extendable delivery jet (see picture page 42)			

Please note that Biotrate is not supplied with reservoir