

Diaphragm Pumps for Air, Gases and Vapours

INNOVATIVE
TECHNOLOGY
WORLDWIDE



NEUBERGER



Series LABOPORT® N 810 FT.18, N 810.3 FT.18 Pumps

LABOPORT® Chemically-resistant Diaphragm Vacuum Pumps

Technical features:

- 100% oil-free transfer
- Pure transfer, evacuation and compression
- Highly compatible with vapours and condensation
- Chemically-resistant
- Therefore suitable for highly aggressive or corrosive gases and vapours
- Maintenance-free
- Environmentally friendly
- Gastight, leakage rate approx. 6×10^{-3} mbar x l/s, not tested in serial production.

The chemically-resistant series N 810 and N 810.3 diaphragm pumps are single- and double-head, dry-running devices used in a wide range of laboratory applications. They transfer and pump down without contamination.

The heart of these very compact pumps is a KNF structured diaphragm. This patented diaphragm was stress-optimized using the Finite Elements method. As a result, we were able to make the pumps smaller while increasing the service life of the diaphragm.

Material in contact with the pumped media

Type/OrderNo.	Pump head	Diaphragm	Valves
N 810 FT.18	PTFE	PTFE-coated	FFPM
N 810.3 FT.18	PTFE	PTFE-coated	FFPM


Technical data:	N 810 FT.18	N 810.3 FT.18
Delivery (l/min) ¹⁾	10	10
Ultimate vacuum (mbar abs.)	100	8
Operating pressure (bar g)	1	1
Connectors for tube (mm)	ID 10	ID 10
Permissible gas and ambient temperature	+5...+40 °C	+5...+40 °C
Voltage/Frequencies	230V/50Hz	230V/50Hz
Mains	IP 44	IP 44
Power P ₁	100 W	90 W
Operating current	0.6 A	0.6 A
Weight	5.9 kg	6.9 kg
Dimensions		
LxHxW (mm)	256/187/146	281/187/140
With thermal switch and power fuse		

Motors with other voltages and frequencies on request.


¹⁾ at atm. pressure

Dimensions and performance characteristics

Pump down time for 10 l receiver

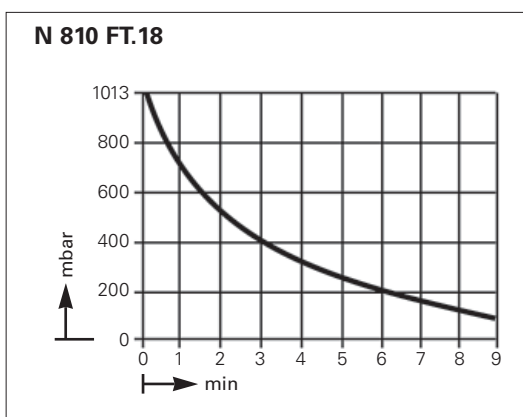


Performance characteristics

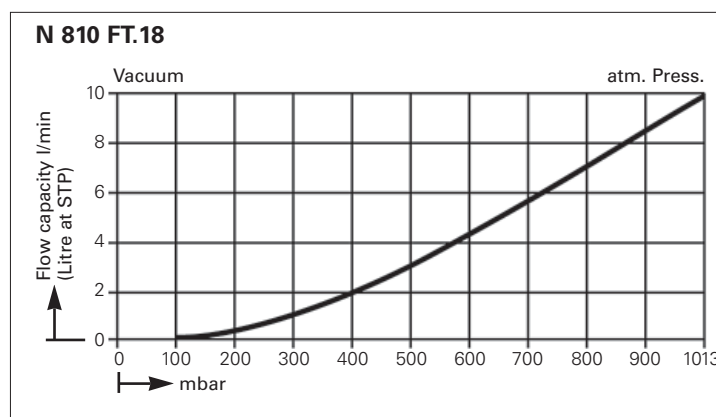


Diaphragm Pumps for Air, Gases and Vapours


Pump down time for 10 l receiver

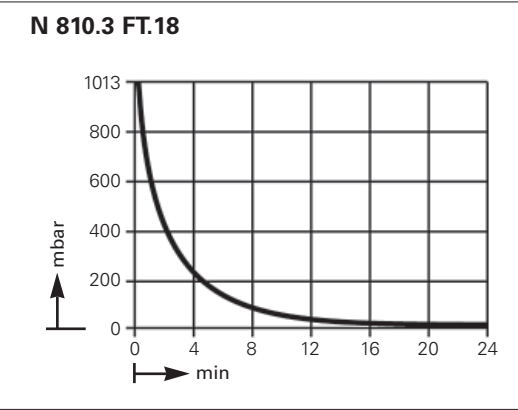

Performance characteristics

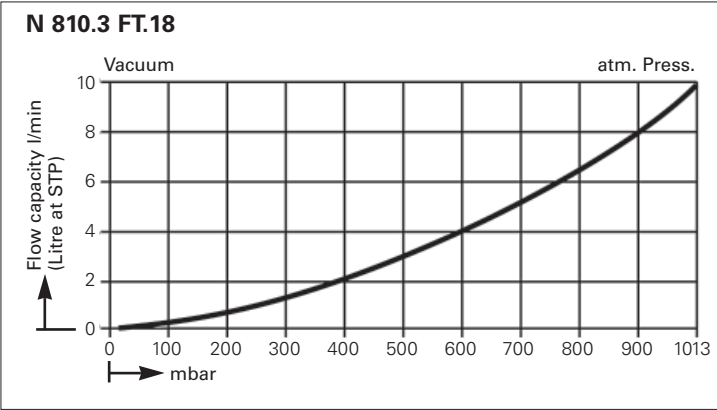

Dimensions (mm)

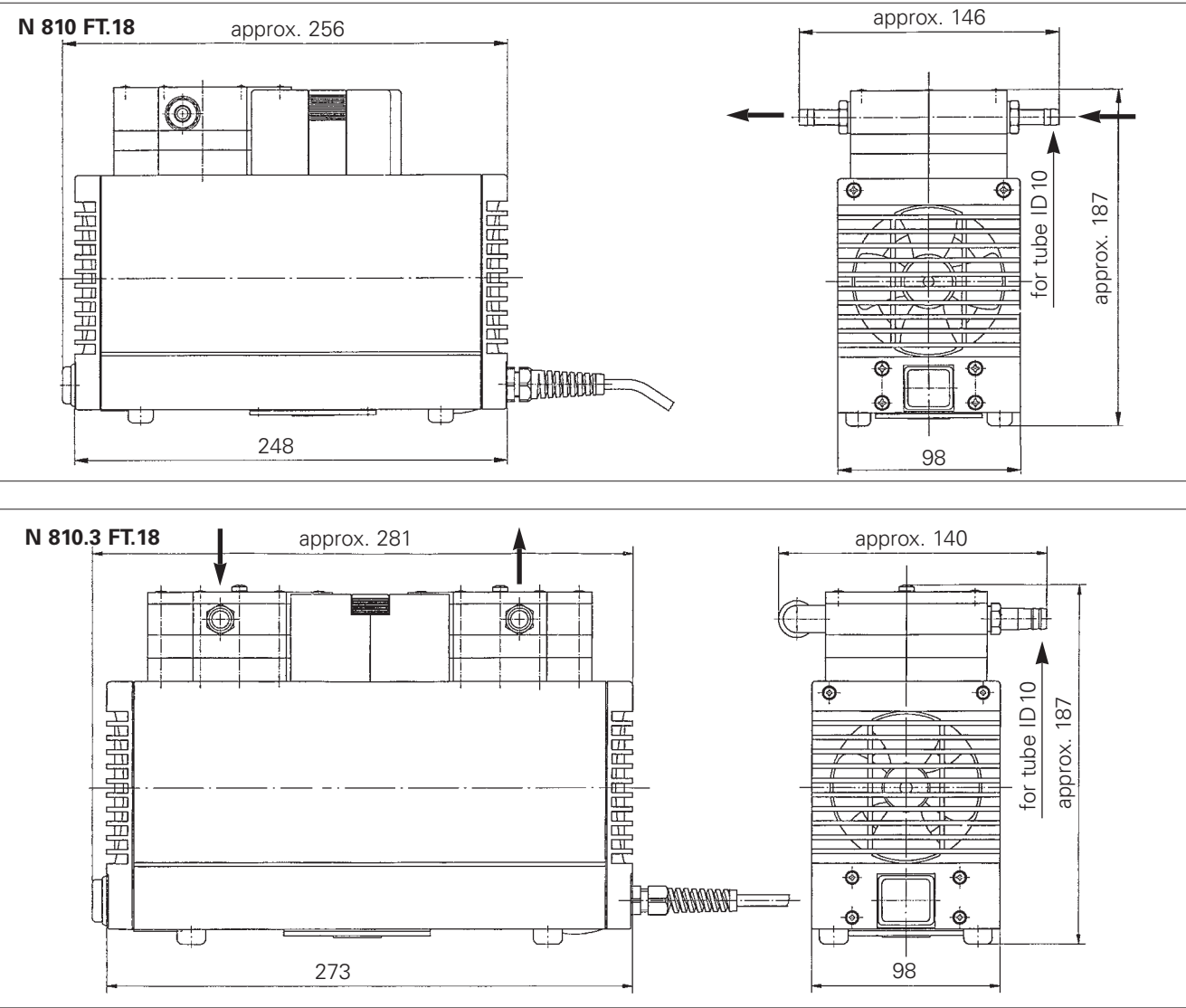

Spare parts

Description	Details	Order No.
Spares kit	for N 810 FT.18	058077
Spares kit	for N 810.3 FT.18	057357



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.