



TECNILABO

HS Code: 84185090

KLAB R700C PL ADV

Laboratory/Pharmacy Refrigerator +5°C

Intended Use Storage of samples/specimens/pharmaceuticals/medications/vaccines

Technical Features

Nominal / Net Storage Volume	700 lt (24,72 cu ft) / 600 lt (21,2 cu ft)
Temperature Range	0°C/+15°C
Temperature Set Point	+5°C
Temperature Uniformity	± 3°C
Temperature Stability	± 1°C
Ambient Operating Temperature	+10°C/+43°C, 50% RH (SN-T Climate Class)
Noise Level	< 42 dBA

Power Supply	220-240V / 50Hz
Standard Electrical Plug	Schuko (electric cable length 3mt)
Power Absorption	2,50 A
Daily Energy Consumption	1,70 kWh/24h (with +20°C ambient temperature)

External Dimensions	W 750 x D 840 x H 2050 mm W 29,52 x D 33,08 x H 80,70 inch.
Internal Dimensions	W 590 x D 670 x H 1510 mm W 23,22 x D 26,38 x H 59,45 inch.
Weight	150 kg

Packing Dimensions	Carton Box + Pallet: W 850 x D 900 x H 2200 mm, 155 kg Wooden Cage: W 1000 x D 1000 x H 2400 mm, 225 kg
---------------------------	--

External Surface	Satin-finished AISI 304 stainless steel
Internal Surface	AISI 304 18/10 stainless steel
External/Internal Edges	Rounded for easy and optimal cleaning
Insulation Type	High-density CFC-free PUR foam
Insulation Thickness	80 mm on sides / 80 mm rear
Door Type	Solid, self-closing system (with opening angle ≤90°), ergonomic handle, key-lock
Door Sealing	Magnetic silicone gasket placed all over the door's frame
Lighting	LED light, with automatic switching on when the door is opened and/or from the display icon


Standard Configuration	- n°3 wire shelves made in PVC-coated steel, adjustable in height (Dimensions W 530 x D 650 mm) with max load 35 kg each - n°4 swivel castors (front castors equipped with brakes) - Porthole (Ø 25 mm) for introducing in the chamber independent probes connected to external monitoring system/BMS - 2,7A backup battery with self-charging circuit (36 hours autonomy in case of power failure) - Volt-free dry contacts (NC-NO) for remote alarms connection to an independent alarm monitoring system /BMS
Max Shelves / Drawers	7 / 7 (or a mix of them)

Cooling System

Cooling	Forced-air cooling, with automatic shut-off of the evaporator fan during door opening
Compressor	Fully sealed hermetic compressor with thermal protector
Refrigerant (CFC/HCFC free)	R290 (gr. 72)
Defrost System	Fully automatic (hot gas) with automatic evaporation of condensation water

PL ADV Control System


Controller	Microprocessor Controller
Display	LED display with 5-keys capacitive keyboard, with 0.1°C resolution
Temperature Probes	n°2 NTC probes (for evaporator / control) + n°1 Pt1000 (for monitoring)
Access Control	3-levels password security system (User/Service/Admin)
Connectivity	In-built Wi-Fi module for remote connection to mobile APP/Cloud
Features	- The main screen displays the temperature set point, functional icons, current time, alarm events, LED light ON/OFF switch - Data logger function with automatic recording of temperature data (at 30-seconds intervals) and events - Flash memory capable of recording data for a minimum period of 1 year - Possibility to mute the alarm notifications - Safety thermostat with anti-freezing (factory pre-set limit at +2°C) and anti-heating (factory pre-set limit at +8°C)
Remote Management and Monitoring with IOS/Android ConneX APP and/or KW Cloud	- Possibility of remotely setting and modifying the device parameters and to view the temperature graph in real time - Possibility to download data in CSV/PDF format via APP on smartphone/tablet/PC (no need for USB) - Automatic sending of push notifications (via APP) and/or email (via Cloud) in the event of an alarm - Possibility of setting automatic sending of weekly temperature and event reports via email
Main Alarms List (Audio/visual notification)	- Power failure - Door open (the alarm interval can be freely set by the user) - High/Low temperature (alarm limits can be set freely by the user) - Temperature probes failure - Back-up battery failure / low charge

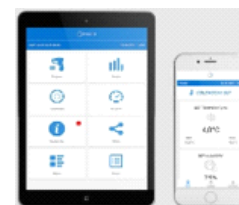



TECNILABO

Standard Features

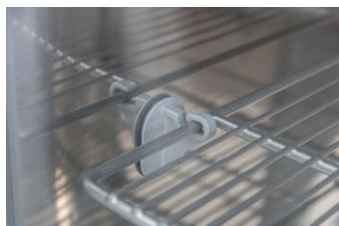

Adjustable wire shelf with fixing system



LED light



Swivel castors with brake



Int/Ext. Porthole



Door switch



Remote alarm contacts



PL ADV controller



Door gasket



Back-up battery

Optional Accessories Available On Request

Additional height-adjustable wire shelf made in PVC-coated steel (n°4 fixing hooks included)
(Shelf's size W 530 x D 650 mm)



Weekly cycle chart disk recorder (including n°52 spare paper disks + n°1 pen)




Perforated drawer (sliding on telescopic guides) made in AISI 304 18/10 stainless steel with plexiglass front panel
(Drawer's size W 535 x D 635 x H 130 mm)



Voltage stabilizer



Perforated drawer (sliding on telescopic guides) made in AISI 304 18/10 stainless steel, lockable for narcotics storage
(Drawer's size W 535 x D 635 x H 130 mm)

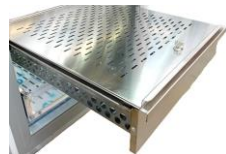

n°2 front leveling feet kit



Plexiglass dividers kit for dividing the drawer into 44 separate sections



Inner power socket (max 2A) with external switch



Pt100 temperature probe enclosed in vial filled with glycol solution (for product simulation)



Packing in wooden cage




Cloud & APP Remote Control



ConneX

FLEX CONTROLLER
MANAGEMENT

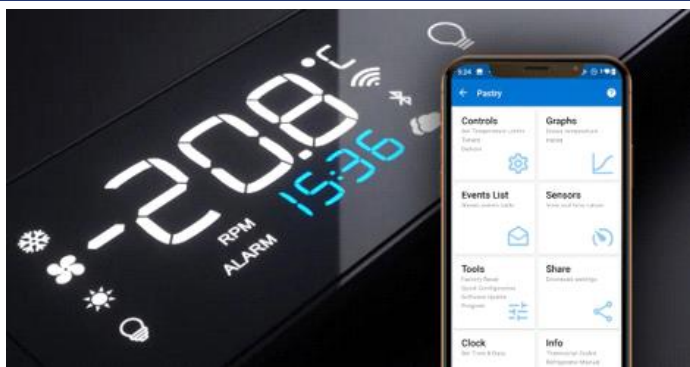

ConneX is the free application (IOS/ANDROID) that allows end users to wirelessly monitor and control PL ADV/HTS controllers installed in refrigeration units. It allows complete control in real time and also remotely of one or more devices, to consult the weekly temperature graph (downloading the data in PDF/Excel format) and to receive alerts via push messages and email.



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.