

Technical Features

Net Storage Volume	708 lt (25,00 cu ft)
2" Cryoboxes Max Capacity	500
1,5/2ml Vials Max Capacity	50.000
Racks Max Capacity	20
Temperature Setting Range	-40°C/-90°C
Temperature Control Range	-40°C/-86°C
Factory T Set Point	-80°C
Ambient Operating Temperature	+10°C/+35°C, 50% RH (N-ST Climate Class)
Noise Level	< 52 dBA

Power Supply	230V ±10% / 50Hz
Standard Electrical Plug	Industrial plug CEE 16A, 2P+E
Max Power Absorption	6A
Daily Energy Consumption*	11,80 kWh/24h




External Dimensions**	W 1060 x D 995 x H 1996 mm W 42,52 x D 39,17 x H 78,58 inch.
Internal Dimensions	W 802 x D 691 x H 1280 mm W 31,57 x D 23,97 x H 50,39 inch.
Weight	340 kg

Packing Dimensions	Carton Box + Pallet: W 1200 x D 1200 x H 2250 mm, 350 kg Wooden Crate: W 1300 x D 1300 x H 2250 mm, 440 kg
---------------------------	---

External Surface	Anti-corrosion hot-dip galvanized steel sheet, covered with PVC film
Internal Surface	AISI 304 18/10 stainless steel
External/Internal Edges	Rounded to facilitate cleaning
Insulation Type	High-density CFC-free PUR foam
Insulation Thickness	140 mm on sides / 140 mm rear / 110 mm outer door
Door Type	Solid, with n°3 heavy-duty side hinges for the max stability
Door Sealing	Triple-layer silicone gasket, heated by the refrigerant to avoid the formation of ice
Handle	Ergonomic with integrated key-lock system (n°2 keys included)

Standard Configuration	n°3 shelves in stainless steel (AISI 304), adjustable in height (adjustment steps of 1.5 cm) with max load of 55kg each n°4 insulated inner doors (inner chamber divided into 4 compartments); inner doors can be removed for cleaning n°4 swivel castors, each one equipped with a leveling foot for a total ground stabilization Porthole (Ø 25 mm) for introducing in the chamber independent probes connected to external monitoring system/BMS External vacuum relief valve (automatic & manual pressure release); the valve can be removed for cleaning operations Back-up battery with self-charging circuit to ensure recording of T and alarm data in case of power failure (48 h autonomy) Volt-free contact for connection to an external alarm monitoring system and BMS
-------------------------------	--

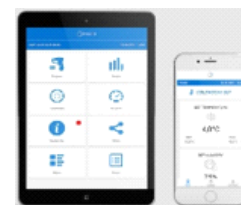

Cooling System

Cooling	Static cooling, with 2-stage cascade refrigeration system; air-cooled condenser with filterless design
Compressor	n°2 hermetic compressors (with thermal protector)
HC Refrigerant (CFC/HCFC free)	1st stage: R1270 (gr. 150) 2nd stage: R170 (gr. 100)
Defrost System	Manual

PL ADV Control System

Controller	Microprocessor Controller
Display	LED display with 5-keys capacitive keyboard, with 0.1°C resolution
Temperature Probes	n°2 NTC probes (for evaporator / temperature control) + n°1 Pt1000 (for monitoring)
Access Control	3-levels password security system (User/Service/Admin)
Connectivity	In-built Wi-Fi module for remote connection to mobile APP ConneX (for IOS & Android)/Cloud The main screen displays the temperature set point, functional icons, current time, alarm events Data logger function with automatic recording of temperature data (at 30-seconds intervals) and events Flash memory capable of recording data for a minimum period of 1 year In-built power stabilizer with linear voltage regulation (± 20V)
Features	Possibility of remotely setting and modifying the device parameters and to view the temperature graph in real time Possibility to download data in CSV/PDF format via APP on smartphone/tablet/PC (no need for USB) Automatic sending of push notifications (via APP) and/or email (via Cloud) in the event of an alarm Possibility of setting automatic sending of weekly temperature and event reports via email
Remote Management and Monitoring with IOS/Android ConneX APP and/or KW Cloud	
Main Alarms List (Audio/visual notification)	- Power failure - Temperature probes failure - Door open (interval can be freely set by the user) - Back-up battery failure / low charge - High/Low T (limits can be set freely by the user) - Dirty condenser / High ambient temperature

Accessories Available On Request

Frontal/Lateral extraction racks made in AISI 304 stainless steel for cryoboxes storage

Additional stainless steel (AISI 304) shelf

Additional internal/external porthole (Ø 25 mm)

Voltage stabilizer

Weekly cycle chart disk recorder (including n°52 spare paper disks + n°1 pen)

CO2 backup system with 12VDC (complete with independent thermal sensor/T controller, 24Ah battery) for activation in case of mechanical/power failure

LN2 backup system with 12VDC (complete with independent thermal sensor/T controller, 24Ah battery) for activation in case of mechanical/power failure

Electrical plug different from standard

Packing in wooden crate



Ergonomic handle with key-lock system



Shelves adjustment system




Insulated inner door




Heated triple-layer silicone gasket



PL ADV controller

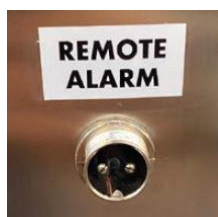

Remote alarm contacts

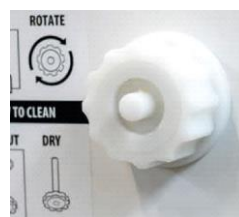

Vacuum relief valve



CO2/LN2 back-up systems



Frontal extraction racks



Lateral extraction racks




* with +20°C ambient temperature

** Exterior dimensions of main cabinet only, excluding handle and other external projections. For complete details, see the dimensional drawings below



TECNILABO

Dimensional Diagram

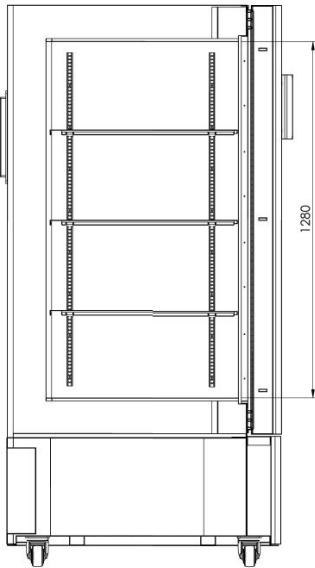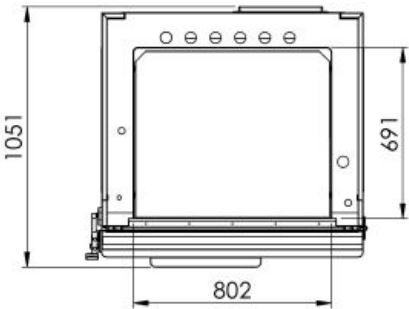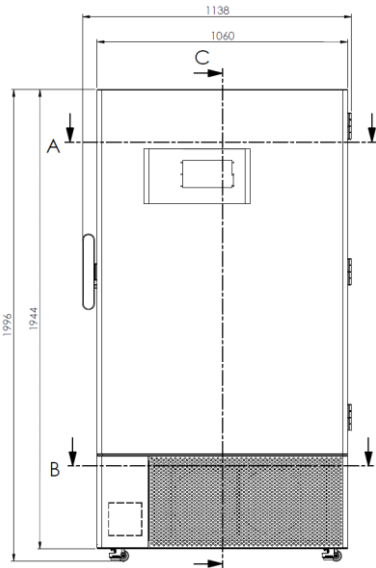



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.