



SL-CLASSIC-UB ◀ 140 cm ▶



Cabinet for acids and alkalis

SL-CLASSIC-UB model SL.060.140.UB.3T

in Light grey (RAL 7035), interior equipment with 4 x pull-out shelf (panel/PP)

Order No. 30616-001-30618



Cabinet for acids and alkalis

SL-CLASSIC-UB model SL.060.059.UB.T

in Light grey (RAL 7035), interior equipment with 2 x pull-out shelf (panel/PP)

Order No. 30622-001-30625

Safe and approved storage of aggressive, non flammable hazardous materials in working areas

Function / construction:

- **robust construction and longevity:** made out of highly resistant melamin resin special plates
- **user-friendly:** access to all containers, 135° opening of the wing door, interior of the cabinet completely visible
- **no unauthorised use:** doors lockable with cylinder lock and locking state indicator (red/green)
- **mobile:** castors with plinth (optionally available), easy moving of the cabinet
- **ventilation:** integrated air ducts ready for connection (DN 50) to a technical exhaust system

Available equipment

- pull-out shelves with removable plastic sump
- **model SL.060.140.UB:** right compartment with integral shelf for material storage



Cabinet for acids and alkalis

SL-CLASSIC-UB model SL.060.110.UB.2T

in Light grey (RAL 7035), interior equipment with 4 x pull-out shelf (panel/PP)


Order No. 30619-001-30621

Technical data


| | | SL.060.140.UB.3T | SL.060.110.UB.2T | SL.060.059.UB.T |
|--|-------|---|-----------------------------------|-----------------|
| External dimensions W x D x H | mm | 1400 x 570 x 600 | 1100 x 570 x 600 | 590 x 570 x 600 |
| Internal dimensions W x D x H | mm | left: 489 x 531 x 562 central: 489 x 531 x 562 right: 281 x 531 x 562 | all compartments: 489 x 531 x 562 | 488 x 531 x 562 |
| Weight without interior equip. | kg | 80 | 60 | 40 |
| Extraction air | DN | 50 | 50 | 50 |
| Air change 30 times | m³/h | 9 | 9 | 4 |
| Diff. pressure (air change 30 times) | Pa | 10 | 10 | 5 |
| Distributed load | kg/m² | 82 | 93 | 151 |
| Maximum load | kg | 100 | 100 | 100 |




SL.060.140.UB.3T Front view

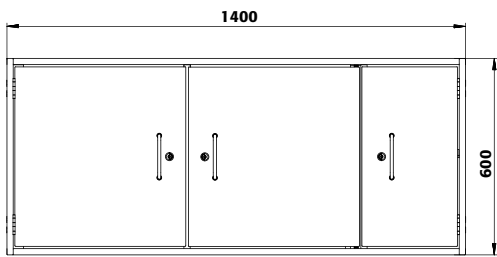


Side view

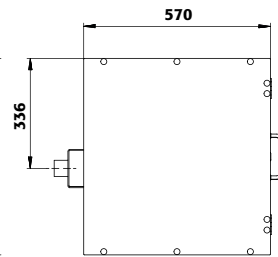


Top view

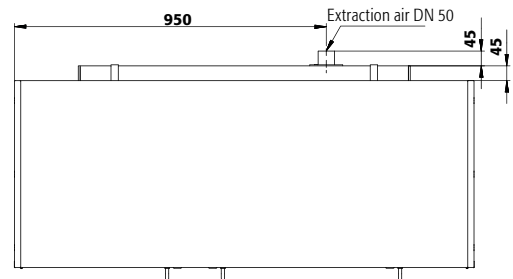


SL.060.110.UB.2T Front view

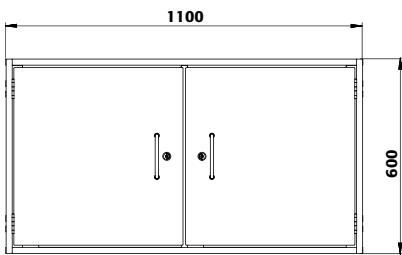


Side view

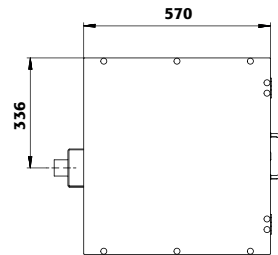


Top view

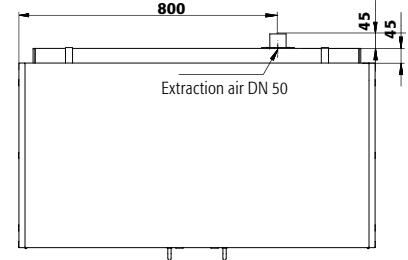


SL.060.059.UB.T Front view

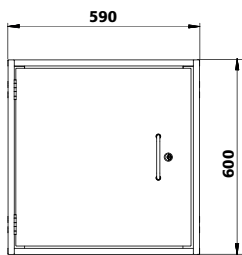

Side view

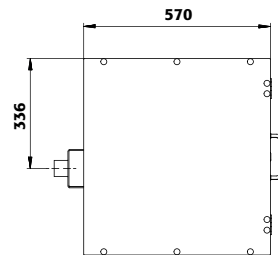

Top view

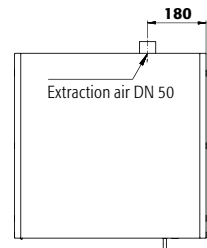

SL-LINE

| Model | Information / equipment | Order No. |
|------------------|--|-----------|
| SL.060.140.UB.3T | without interior equipment, Underbench cabinet, with 3 wing doors | 30616 |
| SL.060.110.UB.2T | without interior equipment, Underbench cabinet, with 2 wing doors | 30619 |
| SL.060.059.UB.T | without interior equipment, Underbench cabinet, left hinged, with wing door | 30622 |
| SL.060.059.UB.TR | without interior equipment, Underbench cabinet, right hinged, with wing door | 30623 |

| Body colour | Door colour | Order No. |
|---------------------|---------------------|--------------|
| Light grey RAL 7035 | Light grey RAL 7035 | 001 4 |

| Interior equipment packages (can only be ordered in combination with the cabinet) | Material | Order No. |
|---|----------|----------------|
| Width 1400 mm | | |
| left: 1x pull-out shelf (V=11.0L) | panel/PP | 30617 2 |
| right: 1x pull-out shelf (V=11.0L) | | |
| left: 2x pull-out shelf (V=11.0L) | panel/PP | 30618 2 |
| right: 2x pull-out shelf (V=11.0L) | | |

| Interior equipment packages (can only be ordered in combination with the cabinet) | Material | Order No. |
|---|----------|----------------|
| Width 1100 mm | | |
| left: 1x pull-out shelf (V=11.0L) | panel/PP | 30620 2 |
| right: 1x pull-out shelf (V=11.0L) | | |
| left: 2x pull-out shelf (V=11.0L) | panel/PP | 30621 2 |
| right: 2x pull-out shelf (V=11.0L) | | |

| Interior equipment packages (can only be ordered in combination with the cabinet) | Material | Order No. |
|---|----------|----------------|
| Width 590 mm | | |
| 1x pull-out shelf (V=11.0L) | panel/PP | 30624 2 |
| 2x pull-out shelf (V=11.0L) | panel/PP | 30625 2 |

| Interior equipment | Material | For cabinet width | 1) | 2) | 3) | Order No. |
|-------------------------------------|------------------------------------|--------------------|-------|-------|-------|----------------|
| pull-out shelf, with PP sump | panel/PP grey | 590, 1100 and 1400 | 25.00 | 18.00 | 16.36 | 30770 2 |
| castors with plinth, height = 30 mm | sheet steel powder-coated RAL 7035 | 1400 | | | | 30374 2 |
| castors with plinth, height = 30 mm | sheet steel powder-coated RAL 7035 | 1100 | | | | 30362 2 |
| castors with plinth, height = 30 mm | sheet steel powder-coated RAL 7035 | 590 | | | | 30478 2 |

| Accessories | Order No. |
|--|----------------|
| extraction module without exhaust air monitoring for wall mounting | 17178 2 |

* usable storage space W x D [mm]

1) load capacity [kg], 2) retention capacity [ltr], 3) maximum container volume in the cabinet [ltr]



Wolflabs

Wolf Laboratories Limited

www.wolflabs.co.uk

Tel: 01759 301142

Fax: 01759 301143

sales@wolflabs.co.uk



Use the above details to contact us if this literature doesn't answer all your questions.

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

