


# Cole-Parmer®

## Cole-Parmer® Tube Rollers, RS-200 and RS-200D

- Multiple rollers saves bench top space and meets standard lab throughput needs
- Featuring removable side panels
- Featuring detachable roller tubes enabling large volume/diameter bottles to fit
- Choice of analogue or digital control
- Greater capacity, ideal for busy large throughput labs, and time saving
- Choice of fixed or variable speed

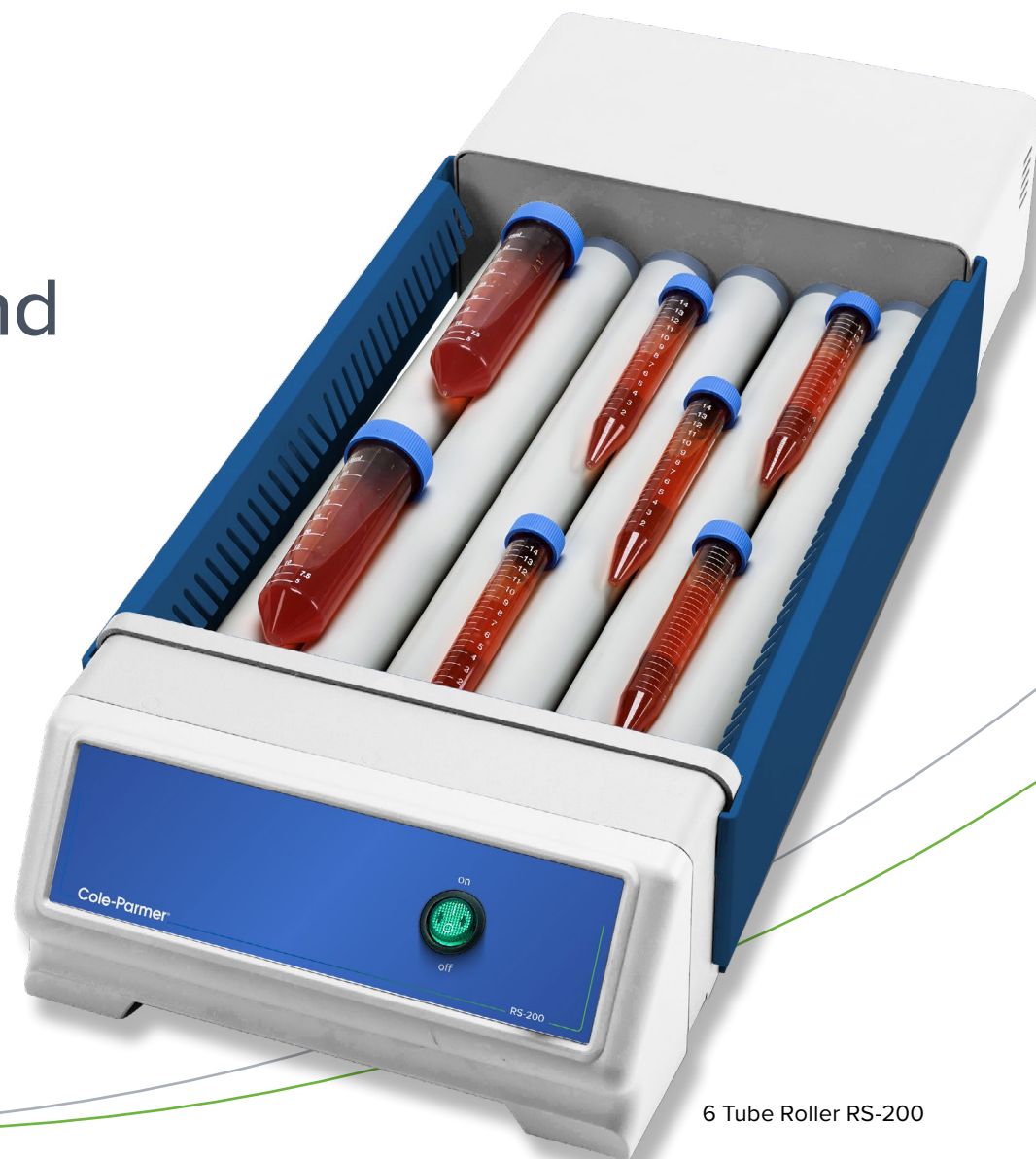

6 Tube Roller RS-200

# Cole-Parmer® Tube Rollers, RS-200 and RS-200D

Cole-Parmer® roller mixers provide a gentle but efficient mixing ideal for mixing blood samples, viscous substances and liquid-solid suspensions where minimum aeration is required. Cole-Parmer® is committed to providing the highest levels of usability and as such has added additional features to the roller mixers. RS-200 roller mixers now feature easy to remove side panels. The side panels provide the user with complete peace of mind that their sample will be mixed effectively during continued use. The side panels ensure the sample bottles do not slide off the mixer during use and therefore provide the highest levels of both safety and functionality.

The roller tubes on Cole-Parmer® roller mixers conveniently detach to accommodate larger vessels and therefore allow the user to mix a wide variety of differently sized tubes and bottles. Cole-Parmer® roller mixers provide the highest levels of versatility as any of the roller tubes can be removed meaning the user can choose the arrangement depending on their needs.



6 Tube Roller RS-200D

## Technical Specification

Specification	RS-200	RS-200D
Number of rollers	6/9	6/9
Speed	33rpm	5 to 60rpm
Amplitude	16mm	16mm
Maximum load	3.0kg	3.0kg
Controls	Analogue	Digital
Timer	No	Yes
Roller size (length x diameter)	340x30mm	340x30mm
Overall dimensions (w x d x h)	595x255x110mm 565x360x110mm	595x255x110mm 565x360x110mm
Net weight	5.1kg & 6.9kg	5.1kg & 6.9kg
Electrical Supply	230V, 50Hz, 50W 120V, 60Hz, 50W	230V, 50Hz, 50W 120V, 60Hz, 50W

## Ordering Information

Description	Ordering Number	Series No.	Model No.	Legacy Sku.
Roller mixer, 6 rollers, analogue control, fixed speed	04653-02	RS-200	RS-200	SRT6
Roller mixer, 6 rollers, analogue control, fixed speed, 120V, 60Hz	04653-01	RS-200	RS-200-120	SRT6/120V/60
Roller mixer, 6 rollers, digital control, variable speed	04653-06	RS-200	RS-200D	SRT6D
Roller mixer, 6 rollers, digital control, variable speed, 120V, 60Hz	04653-05	RS-200	RS-200D-120	SRT6D/120V/60
Roller mixer, 9 rollers, analogue control, fixed speed	04653-04	RS-200	RS-200-L	SRT9
Roller mixer, 9 rollers, analogue control, fixed speed, 120V, 60Hz	04653-03	RS-200	RS-200-L-120	SRT9/120V/60
Roller mixer, 9 rollers, digital control, variable speed	04653-08	RS-200	RS-200D-L	SRT9D
Roller mixer, 9 rollers, digital control, variable speed, 120V, 60Hz	04653-07	RS-200	RS-200D-L-120	SRT9D/120V/60



# Wolf Laboratories Limited

[www.wolflabs.co.uk](http://www.wolflabs.co.uk)

Tel: 01759 301142

Fax: 01759 301143

[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)



**Use the above details to contact us if this literature doesn't answer all your questions.**

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

