


## ISOLATION GLOVE BOX

*for Sorting, Processing, & Opening Mail*

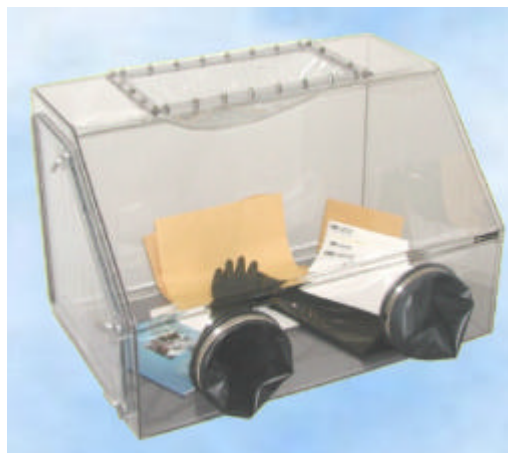

3 ft. Isolation Glove Box for anthrax containment

This Isolation Glove Box is designed to place an added barrier between you and any suspected contaminants, contained within an envelope or package, such as anthrax or small pox. This containment glove box's hermetically sealed interior and neoprene gloves provide the air tight atmosphere you desire when opening tainted mail and the neoprene gloves provide the greatest dexterity available while keeping your skin safe from contamination. The 18" x 18" side door of the glove box allows for larger packages to be passed into the glove box before being opened (side door picture). Finally the patented COY Diaphragm Top allows the user to comfortably enter the chamber through the glove ports without fighting the gloves, and without expelling the interior atmosphere to the room air through a bubbler or pressure relief valve. The COY Isolation Glove Box for anthrax containment is the peace of mind you and your employees have been looking for.

### Why A Diaphragm Top?

As the user enters the gloves their hands are changing the volume inside the glove box. The Diaphragm Top expands and contracts with the changing volume, maintaining a comfortable and uniform internal pressure.



The 18" x 18" Side Door allows larger packages to fit inside the glove box.

## Isolation Glove Boxes Available

Part #	Description	Dimensions
8302-025	2 ft. Glove Box*	24" L x 24" D x 20" H
8302-050	3 ft. Glove Box	36" L x 24" D x 24" H
8302-075	4 ft. Glove Box	48" L x 24" D x 24" H
8302-100	5 ft. Glove Box	60" L x 24" D x 24" H



# Wolf Laboratories Limited

[www.wolflabs.co.uk](http://www.wolflabs.co.uk)

Tel: 01759 301142

Fax: 01759 301143

[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)



**Use the above details to contact us if this literature doesn't answer all your questions.**

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

