

-86°C ULT Freezer

DW-86L51J



Scope of Application :

This product is suitable for pharmaceutical users, the ULT freezer is for the storage of clinical samples and long-term storage of samples in blood stations, hospitals, epidemic prevention stations, scientific research institutes and biological products research institutes as key examples.

Innovative Design

- Ergonomic & User Friendly Design
- Hydrocarbon Energy Saving
- Safe, Secure and Reliable
- Microprocessor Control System
- Tool-free Removable Filter Design
- Multiple Data Ports

Product Advantages




Auto-cascade Hydrocarbon (HC) Refrigeration Technology

Utilizing imported auto-cascade hydrocarbon compressor, with superior refrigeration effect and energy saving



Microprocessor Control System

Microprocessor control, with temperature display precision at 0.1°C. Internal temperature is adjustable from of -40°C ~ -86°C

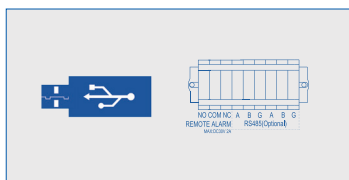

Multiple Ports

Enables users to download historical temperature data for the past 15 years
Equipped with remote alarm function and RS485 port



Filter Design

Tool-free removable filter design, easy and quick to clean the filter



Excellent Doors Seals

Total six gaskets to safeguard the freezer temperature, including a 70mm insulation layer, a 25mm vacuum VIP and 4 layers of silicone seal



Stainless Steel Inner Door

Prevent loss of inside cold air and easy to clean




Ergonomic Design

Unique hook slide, which allows users to open and close with one hand

Internal Lock

Standard internal lock with 4 keys, which can be managed by different people

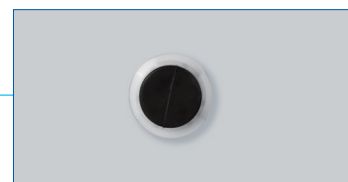

Porthole

One porthole as standard, allows for independent testing of cabinet temperature

Technical Parameters

Model	Power Supply (V/Hz)	Power (W)	Effective Volume (L)	Temperature (°C)	External Dimensions (W*D*H)(mm)	Internal Dimensions (W*D*H)(mm)	Weight (kg)	Sample Size (2ml)
DW-86L51J	220~240/50	320	51	-40~-86	532*640*806	330*481*316	86	3000

*Haier Biomedical reserves the right to change products and specifications without prior notice.



Wolflabs

Wolf Laboratories Limited

www.wolflabs.co.uk

Tel: 01759 301142

Fax: 01759 301143

sales@wolflabs.co.uk



Use the above details to contact us if this literature doesn't answer all your questions.

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

