



OF1202-40MQ3 Oil Free Compressor with Dryer

JUN-AIR®



Features

- 100% oil free
- 100% duty cycle
- Low noise and vibration level
- Easy maintenance
- Long-life, 2-year warranty
- 8000 hours cup life at 8bar
- 4000 hours maintenance interval for the dryer

Description

- Built-in 2-column adsorption dryer
- Pressure dew point down to -40°C
- 1µm air filtration (class 2 for dust according to ISO8573.1:2001)
- Pressure regulator included
- Receiver with internal coating and 5-year warranty


Other Versions

- Without dryer: OF1202-40M
- Other power supply: 400V - 3 phases

Product Specifications


Model	OF1202-40MQ3		
Voltage	V	1x230	1x230
Frequency	Hz	50	60
Motor	HP	2	2
	kW	1.47	1.47
Displacement	l/min	290	328
Air flow @ 8bar	l/min	104	117
	cfm	3.7	4.1
Max. pressure	bar	8	8
	psi	120	120
Max. current	A	8	9
Receiver volume	litres	40	40
	gallon	10.6	10.6
Weight	kg	121	121
	lbs	267	267
Noise level	dB(A)/1m	60	63
Duty cycle		100%	100%

Air Flow / Pressure Curve*

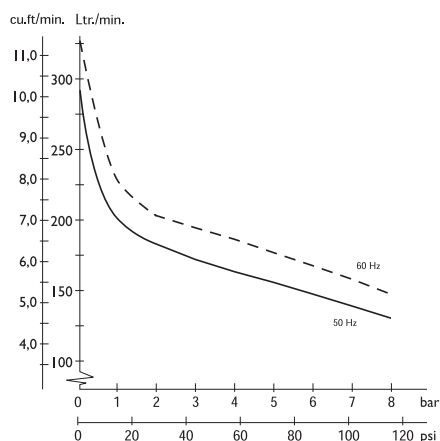

*air flow is reduced by approx. 18-20% on units with dryer

Dimensions*

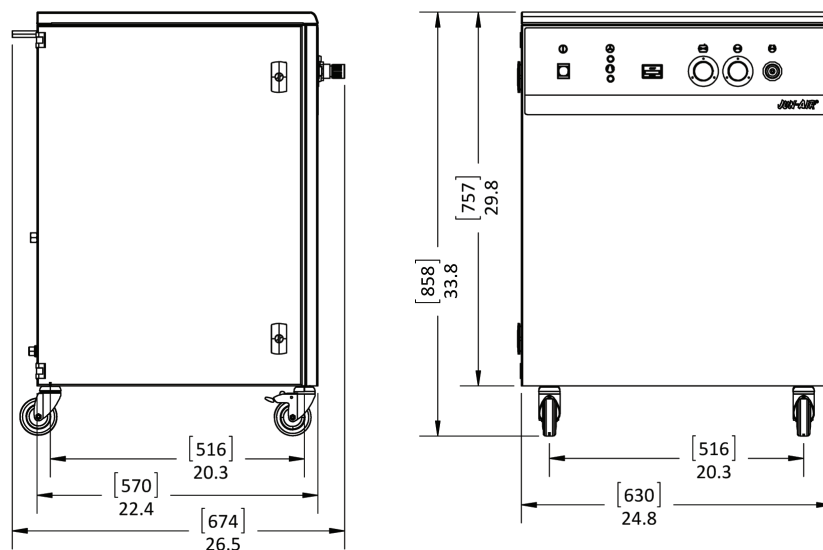

*for unit without accessories or options



Wolf Laboratories Limited

www.wolflabs.co.uk

Tel: 01759 301142

Fax: 01759 301143

sales@wolflabs.co.uk



Use the above details to contact us if this literature doesn't answer all your questions.

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

