

# Diaphragm Pumps for Air, Gases and Vapours

INNOVATIVE  
TECHNOLOGY  
WORLDWIDE



NEUBERGER



## Series LABOPORT® N 840 FT.18, N 840.3 FT.18 Pumps

## LABOPORT® Chemically-resistant Diaphragm Vacuum Pumps

### Technical features:

- 100% oil-free transfer
- Pure transfer, evacuation and compression
- Highly compatible with vapours and condensation
- Chemically-resistant
- Therefore suitable for highly aggressive or corrosive gases and vapours
- Maintenance-free
- Environmentally friendly
- Gastight, leakage rate approx.  $6 \times 10^{-3}$  mbar x l/s, not tested in serial production.

The chemically-resistant series N 840 and N 840.3 diaphragm pumps are single- and double-head, dry-running devices used in a wide range of laboratory applications. They transfer and pump down without contamination.

The heart of these very compact pumps is a KNF structured diaphragm. This patented diaphragm was stress-optimized using the Finite Elements method. As a result, we were able to make the pumps smaller while increasing the service life of the diaphragm.

### Material in contact with the pumped media

Type/OrderNo.	Pump head	Diaphragm	Valves
N 840 FT.18	PTFE	PTFE-coated	FFPM
N 840.3 FT.18	PTFE	PTFE-coated	FFPM


Technical data:	N 840 FT.18	N 840.3 FT.18
Delivery (l/min) <sup>1)</sup>	34	34
Ultimate vacuum (mbar abs.)	100	8
Operating pressure (bar g)	1	1
Connectors for tube (mm)	ID 10	ID 10
Permissible gas and ambient temperature	+5...+40 °C	+5...+40 °C
Mains	230V/50Hz	230V/50Hz
Motor protection	IP 44	IP 44
Power P <sub>1</sub>	180 W	245 W
Operating current	1.5 A	1.5 A
Weight	10.3 kg	12.6 kg
Dimensions LxHxW (mm)	297/226/171	341/226/166
With thermal switch and power fuse		

Motors with other voltages and frequencies on request.


<sup>1)</sup> at atm. pressure

## Dimensions and performance characteristics

### Pump down time for 10 l receiver

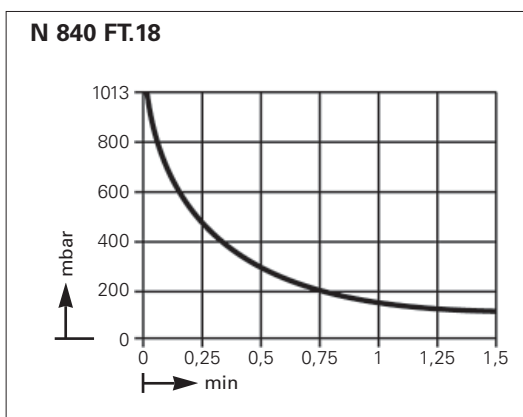


### Performance characteristics

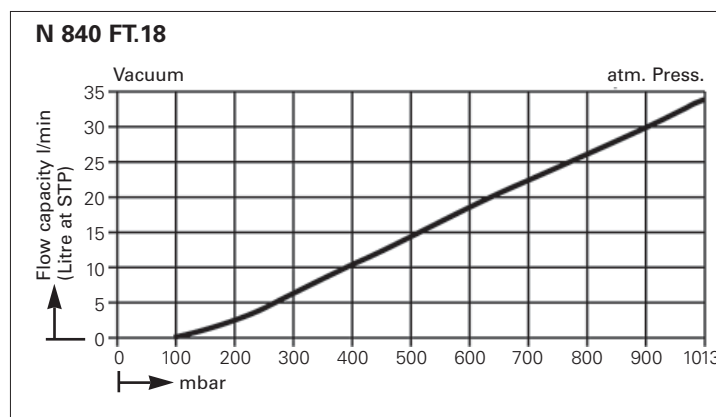


# Diaphragm Pumps for Air, Gases and Vapours


## Pump down time for 10 l receiver

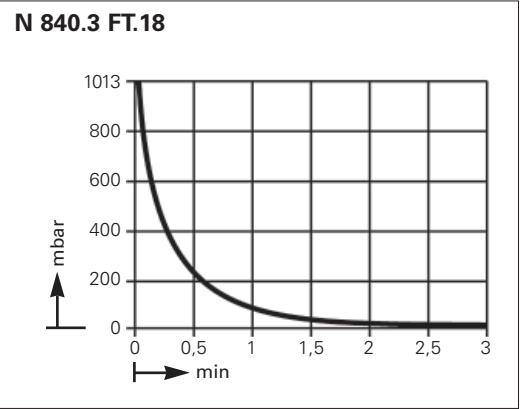

## Performance characteristics

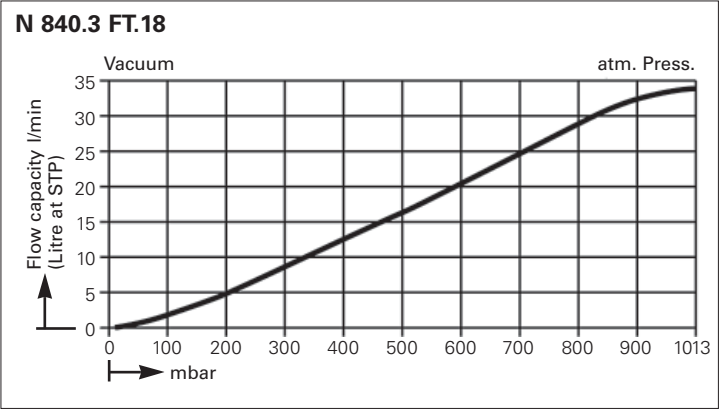

## Dimensions (mm)

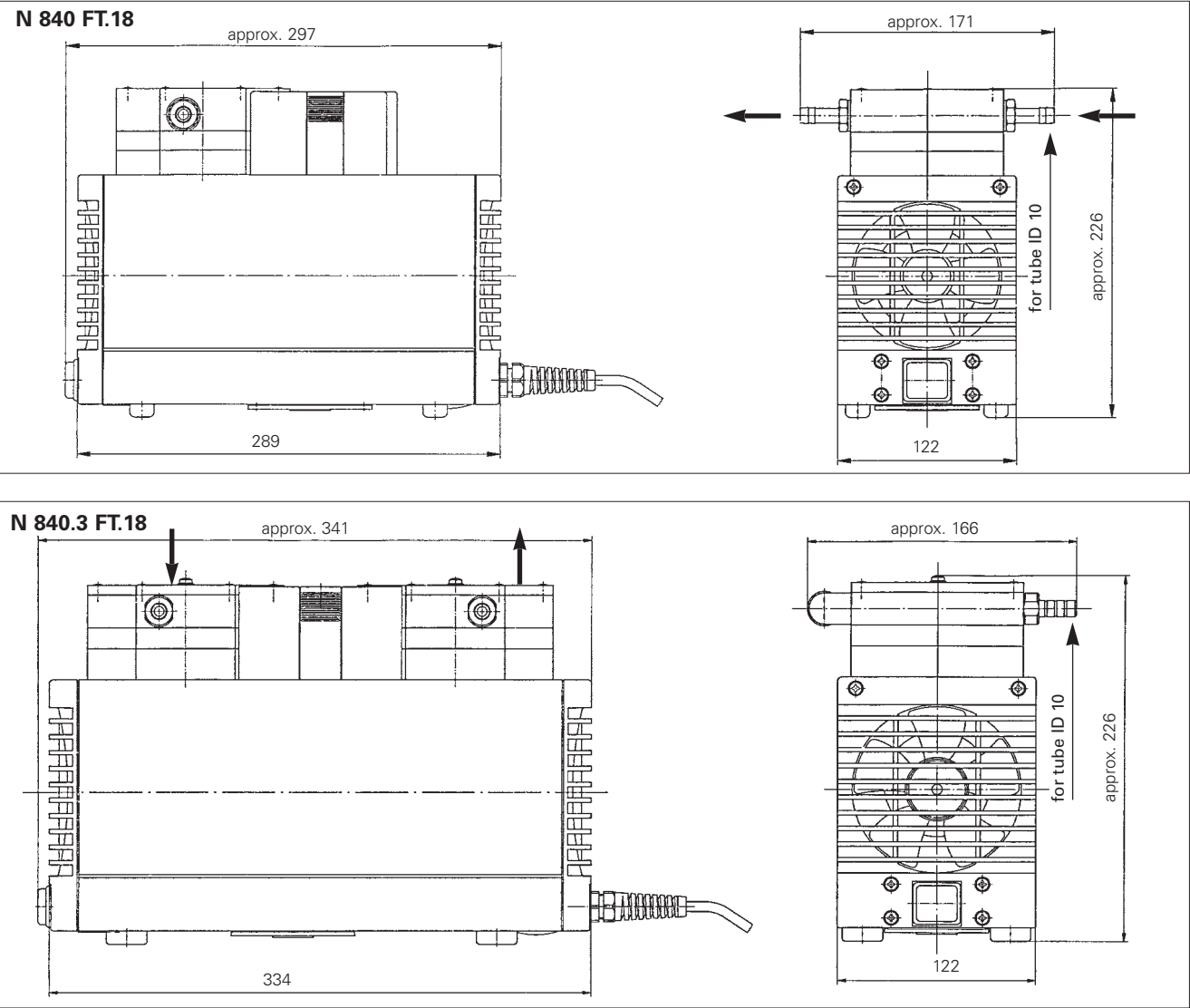

## Spare parts

Description	Details	Order No.
Spares kit	for N 840 FT.18	058079
Spares kit	for N 840.3 FT.18	057359



**Wolflabs**

# Wolf Laboratories Limited

[www.wolflabs.co.uk](http://www.wolflabs.co.uk)

Tel: 01759 301142

Fax: 01759 301143

[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)



**Use the above details to contact us if this literature doesn't answer all your questions.**

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

