


RLVF0217 Bench Top Sparkfree Freezer

Part Number	RLVF0217
Capacity (Litres)	66
Style	Bench Top
Exterior Dimensions (HxWxD)	740 x 450 x 510mm
Interior Dimensions (HxWxD)	590 x 333 x 300mm
Weight	34kg
Energy Consumption	0.19kWh/24hr
Operating Temperature	-18°C to -25°C
Factory Set Point	-20°C
Interior Construction	Powder Coated Steel
Exterior Construction	Powder Coated Steel
Door Style	Self Closing
Door Hinging	Right (can be reversed please ask at time of order)
Shelves	3
Shelf Dimensions (w x d)	2 of 309 x 275mm 1 of 309 x 150mm
Door Lock	✓
Door Open Alarm	✓
Feet	4
Temperature Display	✓
High Temperature Alarm	✓
Low Temperature Alarm	✓
Manual Defrost	✓
Eco Friendly Refrigerant	✓
2 Year Warranty (UK only)	✓

SORRY

No Image

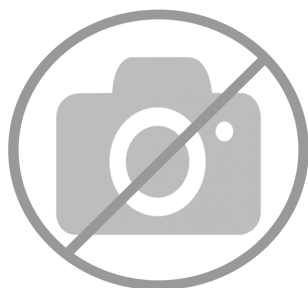

Available

The Labcold Sparkfree range is engineered to store flammable materials safely with no sources of ignition inside the chambers. Anything that could cause a spark is located outside of the chamber.

Warning: Do not store open containers of volatile substances in this freezer.

SORRY

No Image



Available

CLAIM
YOUR

FREE

YEAR

WARRANTY

For extra security all sparkfree models are fitted with a built in lock and digital temperature display. The display features high and low audio and visual temperature alarms as standard in addition to a door open alarm.



Wolflabs

Wolf Laboratories Limited

www.wolflabs.co.uk

Tel: 01759 301142

Fax: 01759 301143

sales@wolflabs.co.uk



Use the above details to contact us if this literature doesn't answer all your questions.

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

