




SORRY

No Image

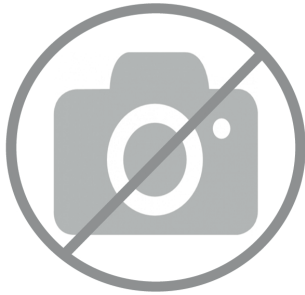

Available

RLVF1517 Sparkfree Freezer

| Part Number | RLVF1517 |
|-----------------------------|---|
| Capacity (Litres) | 406 |
| Style | Free standing |
| Exterior Dimensions (HxWxD) | 1875 x 600 x 700mm |
| Interior Dimensions (HxWxD) | 1690 x 440 x 500mm |
| Net Weight | 97kg |
| Energy Consumption | 3.16kWh/24hr |
| Operating Temperature | -18°C to -25°C |
| Factory Set Point | -20°C |
| Interior Construction | Painted Aluminium |
| Exterior Construction | Painted Steel |
| Door Style | Solid |
| Door Hinging | Right (Can be reversed at time of manufacture) |
| Door Lock | ✓ |
| Door Open Alarm | ✓ |
| Adjustable Shelves | 6 |
| Shelf dimensions (WxD) | 413mm x 438mm |
| Fan Circulated | ✓ |
| Rollers | 4 |
| Temperature Display | ✓ |
| High Temperature Alarm | ✓ |
| Low Temperature Alarm | ✓ |
| Access Port | ✓ |
| Automatic Defrost | ✗ |
| 2 Year Warranty (UK only) | ✓ |

The Labcold Sparkfree range is engineered to store flammable materials safely with no sources of ignition inside the chambers. Anything that could cause a spark is located outside of the chamber.

Warning: Do not store open containers of volatile substances in this freezer.

SORRY

No Image



Available



For extra security all sparkfree models are fitted with a built in lock and digital temperature display. The display features high and low audio and visual temperature alarms as standard in addition to a door open alarm.





Wolf Laboratories Limited

www.wolflabs.co.uk

Wolflabs

Tel: 01759 301142

Fax: 01759 301143

sales@wolflabs.co.uk



Use the above details to contact us if this literature doesn't answer all your questions.

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

