

Shaking Baths



Clifton shaking water baths are ideal for day-to-day use. Available in digital or analogue option. Featuring a broad range of accessory flask clips and racks ensures flexibility for routine and specialised work.

Comprising of a high quality corrosion resistant stainless steel tank, with a bright clean finish housed in a durable painted outer case. Painted surfaces feature anti-bacterial finish with good chemical resistance. Fitted with non slip rubber feet.



Accessory clips and racks can be easily mounted onto the perforated tray which fits onto the reciprocating stainless steel trolley. The trolley can also be easily removed for cleaning.



Analogue Shaking Bath



Analogue shaking water bath features an analogue temperature and shaking speed control. An adjustable over temperature device can be set a few degrees above the set working temperature, protecting samples from damage and over heating.

Digital Shaking Bath



Digital shaking water baths allow temperatures to be easily and rapidly set with the push button, wipe clean control panel, featuring a bold LED display. The high precision PID temperature controller, ensures that accurate results and repeatability is guaranteed.

NE5 Series Shaking Baths



The [Clifton NE5-D shaking water bath](#), features precision PID control, allowing observation of set and actual temperature simultaneously. A broad range of accessory flask clips and racks ensures flexibility for routine and specialized work.

NE5-28 Clifton shaking water bath, features analogue temperature and shaking speed control available with a broad range of accessory trays, flask clips and racks ensures flexibility for routine and specialised incubations.




Comprising of a high quality corrosion resistant stainless steel tank, with a bright clean finish housed in a durable painted outer case. Painted surfaces feature anti-bacterial finish with good chemical resistance.





	NE5 - 10D	NE5 - 28D	NE5 - 28
GENERAL			
Temp. Range	Ambient +5°C - 99°C (-15°C with accessory cooling)		Ambient +5°C - 99°C
Temp. Alarm	Under/over temperature alarm is automatically set 4°C above and 4°C below set point.		Adjustable over temperature device
Temp. Control	Digital touch screen, high precision controller, resolution 0.1°C		Analogue dial pointer
Sensitivity/Uniformity	±0.1°C/±0.1°C @ 45°C (optimal water level 75mm to 120mm)		±0.2°C/+0.1°C @ 45°C (optimal water level 75mm to 120mm)
Heater Power (220-240V)	800W	1000W	
Voltage	230V models 50-60Hz		
Linear Shaking Control	from 0 to 400 Strokes/0 to 200 rpm		
Stroke length	Fixed 20mm Stroke	Variable 40mm Stroke	
Speed Display	Analogue		
Drain (12mm Bore)		Drain Outlet	
DIMENSIONS			
Tank Capacity	10 litres	28 litres	
External (mm)	449w x 332 x 290h	715w x 332d x 290h	
Internal (mm)	242w x 300d x 200h	500w x 298d x 200h	
Perforated Shelf/Trolley (mm)	175w x 270d	376w x 276d	
Water Level Max. (mm)	120		
Water Level Min. (mm)	50		




[Accessories on Page 38 & 39](#)

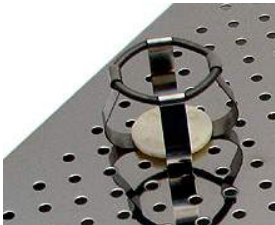

Shaking Baths Accessories

	Stainless Steel Gabled Lid with Insulated Handle	Stainless Steel Flat One Piece lid with Insulated Handle	Clear Polycarbonate Gable Lid with Insulated Handle	Clear Polycarbonate Hinged Gable Lid with Insulated Handle
SHAKING BATHS				
NE5-10D	SL1-8	Not Available	GL1-8AF	GL1-8
NE5-28 & NE5-28D	SL1-22	LD4-8	GL1-22AF	GL1-22

	Stainless Steel Flask Clips FC0025 : 25ml FC0050 : 50ml FC0100 : 100ml FC0250 : 250ml FC0500 : 500ml FC1000 : 1000ml FC2000 : 2000ml	Stainless Steel Perforated Tray	Stainless Steel Test Tube Rack Tray	Drain Tap
SHAKING BATHS				
NE5-10D	Attach selected clips to perforated tray supplied with Bath	SA01773	Not Required - attached test tube rack directly to perforated tray supplied with Bath	BX0953
NE5-28 & NE5-28D	Attach selected clips to perforated tray supplied with Bath	SBF01804	STT004	BX0953

<p>Stainless Steel Flask Clip Tray</p> <p>SFT0025 : 39 x 25ml SFT0050 : 32 x 50ml SFT0100 : 18 x 100ml SFT0250 : 12 x 250ml SFT0500 : 8 x 500ml SFT1000 : 5 x 1000ml SFT2000 : 2 x 2000ml</p>	<p>Stainless Steel Test Tube Racks</p> <p>6870 : 26 holes x 17mm DIA 6871 : 16 holes x 26mm DIA 6872 : 36 holes x 13mm DIA 6873 : 18 holes x 19mm DIA 6900 : 12 holes x 32mm DIA</p>	<p>Polypropylene Spheres: Pack of 200</p>	<p>Universal Spring Tray</p>
			
Not Available	2 Rack From Above Range	BP0368 : 1 Pack	Not Available
1 Clip Tray from Above Range	4 Rack From Above Range	BP0368 : 2 Packs	SFT9000

Universal Accessories

CAL	Calibration Record
BX0616	Drain Syphon
BKF-1	Barkeepers Friend
TC-1	Thermometer clip with bent stem thermometer
DC1-300	Immersion cooler for temperature to -15°C (not for use with NE5-28) more details on page 42
NE5-MB	Mountain Bracket for use with DC1-300



Wolflabs

Wolf Laboratories Limited

www.wolflabs.co.uk

Tel: 01759 301142

Fax: 01759 301143

sales@wolflabs.co.uk



Use the above details to contact us if this literature doesn't answer all your questions.

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

