

SciQuip Incu S-Series


- Capacity 50ltr – 270ltr
- Temp ambient +5 to 65°C
- Microprocessor control
- Fully programmable
- Electro-polished chamber & shelves
- Timer or continuous
- Audible & visual alarms
- Password protection
- Large LCD screen
- Separate overheat safety protection
- Inner glass door



The SciQuip Incu S-series incubators are ideal for microbiological work, pharmaceuticals, cosmetics, food and general testing. The S-series incubators are fan assisted providing gentle air flow and excellent temperature uniformity. All parameters are displayed in the large LCD screen for fast easy operation and read-out. The SciQuip Incubators can store up to 9 programs of 18 steps, each program can be made to repeat from 1 to 99 times automatically. All S-Series incubators stand on feet, with the exception of the Incu-270S, which is on lockable castor wheels.



Control Panel

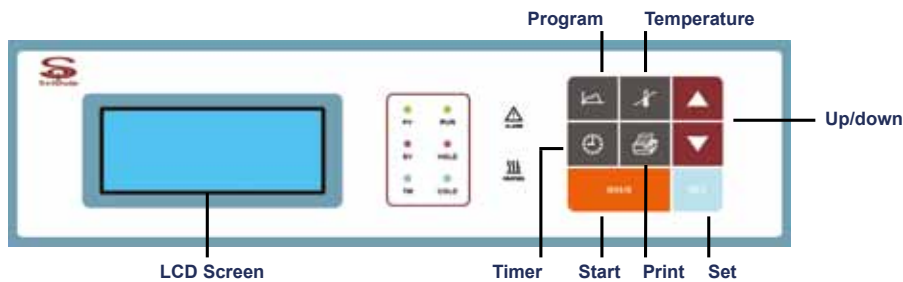

Technical Specifications

Model	Incu-50S	Incu-80S	Incu-120S	Incu-160S	Incu-270S
Capacity (ltrs)	50	80	120	160	270
Temperature range (°C)	ambient +5 to +60	ambient +5 to +60	ambient +5 to +60	ambient +5 to +60	ambient +5 to +60
Temperature set division (°C)	+/-0.1	+/-0.1	+/-0.1	+/-0.1	+/-0.1
Temperature uniformity (°C)	+/- 0.5	+/- 0.5	+/- 0.5	+/- 0.5	+/- 0.5
Timer (mins)	999, continuous	999, continuous	999, continuous	999, continuous	999, continuous
Shelves supplied	2	2	2	2	Max shelves
Max shelves	4	4	4	4	4
Inner dimensions (mm)	350x350x410	400x400x500	450x450x600	500x500x650	600x550x820
Outer dimensions (mm)	470x465x760	520x515x850	570x565x950	620x615x1000	720x665x1170
Packed dimensions (mm)	590x585x930	640x635x1020	690x685x1120	740x735x1170	840x785x1340
Weight (kg)	33	40	51	57	90
Power V/Hz	220-230 50/60	220-230 50/60	220-230 50/60	220-230 50/60	220-230 50/60
Warranty	2yrs	2yrs	2yrs	2yrs	2yrs



Wolflabs

Wolf Laboratories Limited

www.wolflabs.co.uk

Tel: 01759 301142

Fax: 01759 301143

sales@wolflabs.co.uk



Use the above details to contact us if this literature doesn't answer all your questions.

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

