


Key Features

- Compact unit
- Ideal for microtitre plates and petri dishes
- Laminar flow air circulation
- Digital setting and display of temperature
- Convenient upward opening door
- Stainless steel interior



SI19

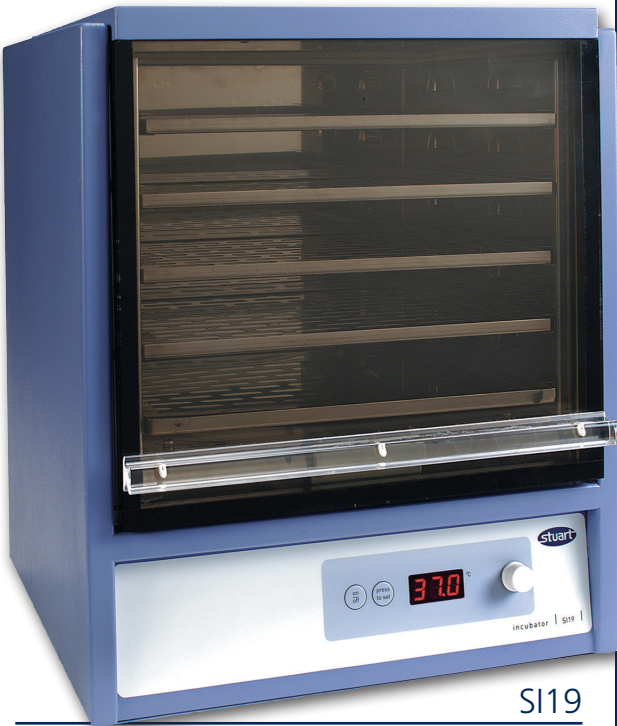

SI19

SI19



Incubator, for microtitre plates,

Specially designed for use with microtitre plates, which can be difficult to incubate. In conventional incubators, when 4 or more plates are placed on a shelf, they block the air circulation, which can seriously impair temperature distribution.

The Stuart® SI19 Incubator overcomes this problem with forced air circulation, which blows preheated air across the shelves rather than just through them. Microtitre plates are also prone to drying out during incubation. The SI19 overcomes this with a water tray to increase the humidity within the chamber.

The bright LED display indicates temperature and the encoder control allows rapid and accurate temperature setting. An integral over-temperature protection system automatically tracks the set temperature and controls the heater in the event of a fault. Error codes are displayed if a fault has occurred or if the set temperature has been interrupted.

An acrylic upwards opening door provides excellent access to the interior. The stainless steel chamber features 4 shelves each accommodating 6 plates (or 12 stacked in 2 layers). Shelves are captive and slide out for easy loading.

Technical Specification

Nominal capacity	20 litres
Temperature range	Ambient + 8°C to 80°C
Temperature fluctuation (at 37°C)	±0.5°C
Temperature variation	±0.5°C
Display type	LED
Display resolution	0.1°C
Internal dimensions, mm (w x d x h)	250 x 230 x 200
Overall dimensions, mm (w x d x h)	380 x 380 x 435
Net weight, kg	24
Electrical supply	230V, 50Hz, 280W
IP Rating	30

Ordering Information

Model	Description
SI19	Incubator, microtitre plate



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.