



Rotavapor® R-300

Convenient and efficient rotary evaporation

The R-300 meets the highest expectations in convenience and versatility. Its modular design allows extension of the Rotavapor® R-300 to a fully integrated system where a central interface regulates each component.

Convenient

Ergonomic handle, easy to use Combi-Clip and intuitive user interaction



Reproducible

Precise operation control of all process parameters



Extendable and versatile

Upgrade your individual Rotavapor® R-300 to a fully integrated system at any time





Rotavapor® R-300

Key features and advantages



reddot award 2016
winner



DESIGN
AWARD
2016



GERMAN
DESIGN
AWARD
SPECIAL
2016



Configure your individual Rotavapor® system:
www.buchi.com/Rotavapor





Ergonomic handle

Ergonomic handle to conveniently change the position of the lift. Electronic and manual version available.



Direct rotation speed control

Immediate adjustment of rotation speed with intuitively positioned control knob in case of an urgent or unexpected event.



Combi-Clip

The multifunctional Combi-Clip with snap lock mechanism allows connection of the evaporating flask with one hand.



Integrated display

The heating bath temperature, rotation speed and lift position are shown on the integrated heating bath display.



Additional interface

Control and regulate all process parameters conveniently with the optional Interface I-300 or I-300 Pro.




Glass assemblies

Seven different glass assemblies enable a wide range of applications, solvents and sample properties.

Heating bath options

Choose between two different heating baths. The B-301 for 1 L flask and 95 °C maximum temperature features a very short heat-up time. For flasks up to 5 L and 220 °C maximum temperature, the B-305 covers a broader range of applications.

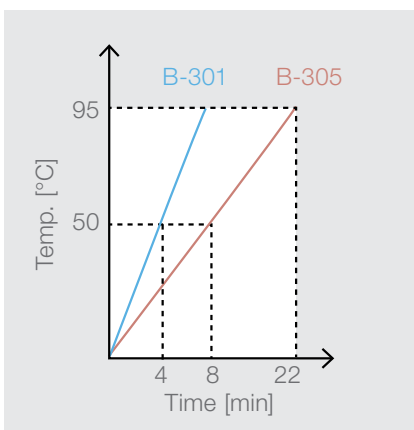

B-301



B-305

4 min to 50 °C

8 min to 50 °C

P_{max} 1100 W

P_{max} 1300 W

1 L flask

5 L flask

T_{max} 95 °C

T_{max} 220 °C



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.