

# NE4 - SB Universal Heated Thermostirrer

*Clifton Range®*




Clifton Range NE4-SB features a digital PID temperature controller which provides even water circulation and temperature stability.

Designed to be fully compatible and easy to fit onto an array of tanks and vessels. Built in the United Kingdom using stainless steel and other premium quality materials the NE4-SB is backed with a three year warranty. Featuring a low-level float device, which ensures heating is switched off if the liquid falls below the recommended level.

## NE4 - SB

GENERAL	
Temp. Range	Ambient + 5°C - 99°C (-20°C with accessory cooling)
Temp. Alarm	Under/over temperature alarm is automatically set 4°C above and 4°C below set point
Temp. Control	PID Digital, LED display, resolution 0.1°C
Sensitivity	±0.05°C
Timer	Visual indicator and audible end message - Variable 00:01 to 99:59 hh:mm
Heater Power (220-240V)	2000W
Voltage	230V models available 50-60Hz
DIMENSIONS	
Max Tank Capacity	60 Litre
External (mm)	122w x 240d x 335h mm
Tank mounting (mm)	Double clamp fitting accommodates tank max wall thickness 40mm
Tank water Levels (mm)	Example : with overall tank depth 150mm / max water level 130mm / min 80mm (see below)

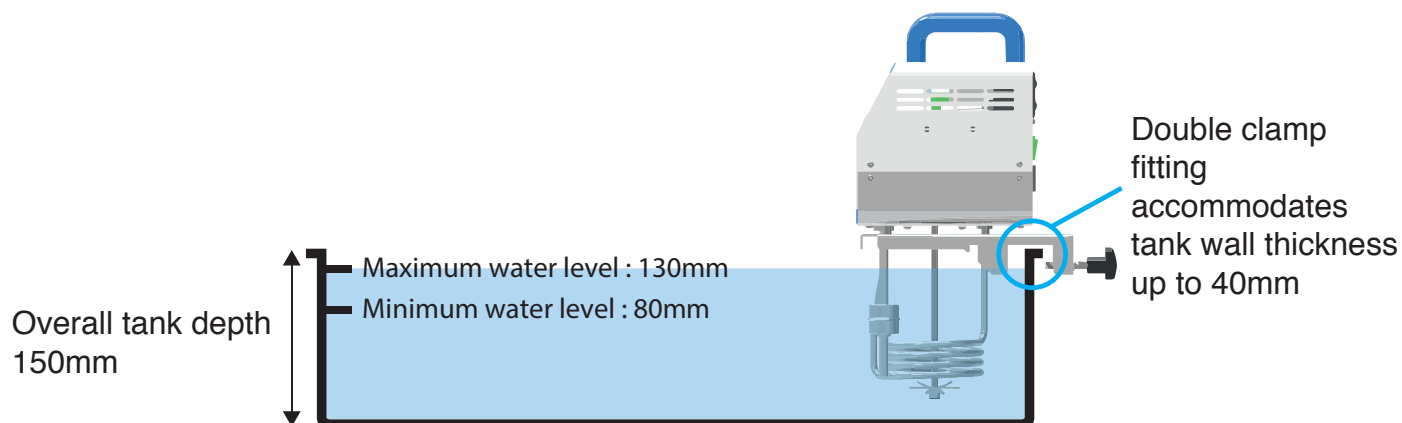



# WolfLabs

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

**[www.wolflabs.co.uk](http://www.wolflabs.co.uk)**

**Tel : 01759 301142**

**Fax : 01759 301143**

**[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)**

Please contact us if this literature doesn't answer all your questions.