

Smart Series (Square Racks)



Scope of Application:

The smart, IoT and cloud management system monitors temperature and liquid levels simultaneously to provide accurate and real-time information on the critical parameters to ensure ultimate sample safety.

Innovative Design

- IoT and Cloud Management
- Multiple Level of Protections
- Low Energy Consumption
- Temperature and Level Display
- Updated Design
- High Performance Stability
- New Lock Design


Key Features

- 5 models from 2,400 to 6,000 cryovial capacity
- 5 years vacuum warranty
- Durable aluminium construction
- Real time monitor of temperature and level
- Cloud data storage for traceability
- Low consumption rate and high performance stability
- New lock design

Product Advantages

Colour identification for Rack handles

Rack handles are colour-coded for ease of separating storage zones and managing samples.

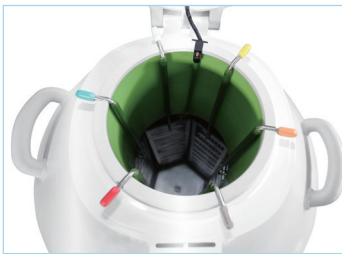


Double lock and double control design

The new double lock and double control design requires two people to open at the same time to ensure sample safety.



Dual monitoring respectively for temperature and liquid level

Storage temperature and liquid level are automatically monitored by a high precision controller. Accurate and real-time information ensures samples safety.

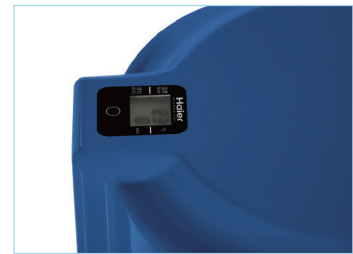

Multilayer protection for ultimate safety

Temperature and level are monitored in real time to safeguard the storage environment. Alarms can be transmitted through email, IM and wechat.


Low consumption rate and high performance stability

Automated manufacture ensures build quality and reliable vacuum which delivers stable and uniform temperature performance as well as low LN₂ consumption rates.



Accessories



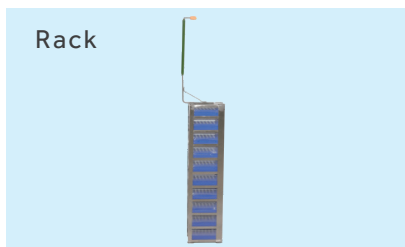


Rack



Roller Base



Technical Parameters

Model	YDS-65-216	YDS-95-216	YDS-115-216	YDS-145-216	YDS-175-216
Picture					
Authentication	CE	CE	CE	CE	CE
Maximum Storage Volume					
No. of canister	6	6	6	6	6
2ml Storage Vials (100/box)	2400	3000	3600	4800	6000
Layers in Rack	4	5	6	8	10
5ml Storage Vials (81/box)	972	972	1458	1944	2430
Layers in Rack	2	2	3	4	5
Function					
Liquid Nitrogen Volume (L)	65	95	115	145	175
Static Evaporation* (L/Day)	0.78	0.94	0.94	0.94	0.95
Static Holding Time** (Day)	83	101	122	154	184
Dimensions					
Neck Opening (mm)	216	216	216	216	216
Overall Height (mm)	760	825	895	995	1110
Overall Diameter (mm)	681	681	681	681	681
Empty Weight (kg)	41	44	48	53	58
Full Weight (kg)	108	121	141	171	199.5
Shipping Weight (kg)	65.5	69.2	75.9	85.2	94.5

*Static Evaporation: The amount of liquid nitrogen lost by natural evaporation within 24 hours after the container is allowed to reach thermal equilibrium under the specified filling rate, which is converted into the evaporation value under the standard state (atmospheric pressure 101.325kPa, ambient temperature 20°C).

**Static Holding Time: under the standard state (atmospheric pressure 101.325kPa, ambient temperature 20°C), the container is filled with liquid nitrogen pre-cooled and reaches thermal equilibrium, and the liquid nitrogen filling rate reaches 100%, the number of days when the liquid nitrogen is completely evaporated on standing.

Product List

Model	Product	Product Model	Description and Accessories
YDS-65-216	Smart Liquid Nitrogen Container (Smart Series)	YDS-65-216	Includes 6 racks(4-tier), cryoboxes(9*9), lockable cap, roller base is optional
YDS-95-216	Smart Liquid Nitrogen Container (Smart Series)	YDS-95-216	Includes 6 racks(5-tier), cryoboxes(9*9), lockable cap, roller base is optional
YDS-115-216	Smart Liquid Nitrogen Container (Smart Series)	YDS-115-216	Includes 6 racks(6-tier), cryoboxes(9*9), lockable cap, roller base is optional
YDS-145-216	Smart Liquid Nitrogen Container (Smart Series)	YDS-145-216	Includes 6 racks(8-tier), cryoboxes(9*9), lockable cap, roller base is optional
YDS-175-216	Smart Liquid Nitrogen Container (Smart Series)	YDS-175-216	Includes 6 racks(10-tier), cryoboxes(9*9), lockable cap, roller base is optional
	Roller Base	YSC-685-5W	Suitable for all 65L-175L aluminium tank, equipped with 5 universal wheels

*Haier Biomedical reserves the right to change products and specifications without prior notice.



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.