


## Compact laboratory balance KERN PCD

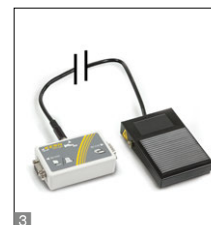

## High-resolution precision balance with removable display for maximum flexibility

### Features

- Laboratory balance with separate platform: Ideal when working in the glove bag or fume cupboard. Particularly practical for weighing toxic, volatile or contaminated substances
- **PRE-TARE function** for manual subtraction of a known container weight, useful for checking fill-levels
- **Freely programmable weighing unit**, e.g. display direct in special units such as length of thread g/m, paper weight g/m<sup>2</sup>, or similar
- **Level indicator and levelling feet** for precise levelling of the scale, fitted as standard, to give the most accurate weighing result
- **1 Draught shield** standard for models with weighing plate size **A**, weighing space W×D×H 146×146×80 mm
- **Protective working cover** included with delivery


### Technical data

- Large backlit LCD display, digit height 21 mm
- Dimensions weighing surface  
**A** Ø 105 mm  
**B** W×D 160×160 mm, see larger picture
- Weighing plate material  
**A** plastic, with conductive lacquer  
**B** stainless steel
- Dimensions of display device W×D×H 140×46×75 mm
- Optional battery operation, 9 V block not included, operating time up to 12 h, AUTO-OFF function to preserve the battery
- Overall dimensions W×D×H 165×280×75 mm
- Cable length of display device approx. 1,2 m
- Net weight approx. 1,2 kg
- Permissible ambient temperature 5 °C/35 °C


### Accessories

- **Protective working cover** over the display device, scope of delivery: 5 items, KERN PCD-A05S05
- **2 Stand** to elevate display device, height of stand approx. 250 mm, KERN PCD-A03
- **Rechargeable battery pack internal**, operating time up to 24 h without backlight, charging time approx. 10 h, KERN PCD-A04
- **3 Foot switch**, ideal when the application requires two free hands. TARE or PRINT function can be selected. Scope of delivery: foot switch, junction box, connection cable. For the PRINT function you will need the RS-232 interface cable, KERN YKF-01
- **RS-232/Ethernet adapter** for connection to an IP-based Ethernet network, KERN YKI-01
- **Individual header data:** the free software can be used to print 4 header lines on the printout when using printers 911-013, YKN-01, YKB-01N, YKE-01 and YKC-01 (in combination with YKI-02)
- Further details, plenty of further accessories and suitable printers see *Accessories*

#### STANDARD

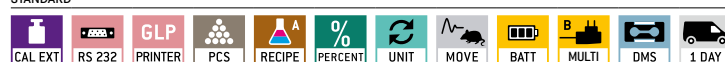



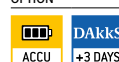














#### OPTION



Model	Weighing capacity [Max] g	Readability [d] g	Reproducibility g	Linearity g	Weighing plate	Quality code	Option	
							DAkkS Calibr. Certificate	
<b>KERN</b>						<b>QUALITY</b>	<b>DAkkS KERN</b>	
PCD 250-3	250	0,001	0,002	± 0,005	<b>A</b>	AA	963-127	
PCD 300-3	350	0,001	0,002	± 0,005	<b>A</b>	AA	963-127	
PCD 2500-2	2500	0,01	0,02	± 0,05	<b>B</b>	BA	963-127	
PCD 3000-2	3500	0,01	0,02	± 0,05	<b>B</b>	BA	963-127	
PCD 6K-4	6000	0,1	0,1	± 0,3	<b>B</b>	BA	963-128	
PCD 10K0.1	10000	0,1	0,1	± 0,3	<b>B</b>	BA	963-128	
PCD 10K-3	10000	1	1	± 3	<b>B</b>	BA	963-128	

## Pictograms

	<b>Internal adjusting:</b> Quick setting up of the balance's accuracy with internal adjusting weight (motordriven)		<b>KERN Communication Protocol (KCP):</b> It is a standardized interface command set for KERN balances and other instruments, which allows retrieving and controlling all relevant parameters and functions of the device. KERN devices featuring KCP are thus easily integrated with computers, industrial controllers and other digital systems		<b>Protection against dust and water splashes IPxx:</b> The type of protection is shown in the pictogram.
	<b>Adjusting program CAL:</b> For quick setting up of the balance's accuracy. External adjusting weight required				<b>Stainless steel:</b> The balance is protected against corrosion
	<b>Easy Touch:</b> Suitable for the connection, data transmission and control through PC, tablet or smartphone				<b>Suspended weighing:</b> Load support with hook on the underside of the balance
	<b>Memory:</b> Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.		<b>GLP/ISO log:</b> The balance displays serial number, user ID, weight, date and time, regardless of a printer connection		<b>Battery operation:</b> Ready for battery operation. The battery type is specified for each device
	<b>Alibi memory:</b> Secure, electronic archiving of weighing results, complying with the 2014/31/EU standard.		<b>GLP/ISO log:</b> With weight, date and time. Only with KERN printers		<b>Rechargeable battery pack:</b> Rechargeable set
	<b>Data interface RS-232:</b> To connect the balance to a printer, PC or network		<b>Piece counting:</b> Reference quantities selectable. Display can be switched from piece to weight		<b>Universal mains adapter:</b> with universal input and optional input socket adapters for A) EU, CH; B) EU, CH, GB, USA; C) EU, CH, GB, USA, AUS
	<b>RS-485 data interface:</b> To connect the balance to a printer, PC or other peripherals. Suitable for data transfer over large distances. Network in bus topology is possible		<b>Recipe level A:</b> The weights of the recipe ingredients can be added together and the total weight of the recipe can be printed out		<b>Mains adapter:</b> 230V/50Hz in standard version for EU. On request GB, USA or AUS version available
	<b>USB data interface:</b> To connect the balance to a printer, PC or other peripherals		<b>Recipe level B:</b> Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display		<b>Power supply:</b> Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, USA or AUS on request
	<b>Bluetooth* data interface:</b> To transfer data from the balance to a printer, PC or other peripherals		<b>Recipe level C:</b> Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display, multiplier function, adjustment of recipe when dosages are exceeded or barcode recognition		<b>Weighing principle: Strain gauges</b> Electrical resistor on an elastic deforming body
	<b>WLAN data interface:</b> To transfer data from the balance to a printer, PC or other peripherals				<b>Weighing principle: Tuning fork:</b> A resonating body is electromagnetically excited, causing it to oscillate
	<b>Control outputs (optocoupler, digital I/O):</b> To connect relays, signal lamps, valves, etc.		<b>Totalising level A:</b> The weights of similar items can be added together and the total can be printed out		<b>Weighing principle: Electromagnetic force compensation</b> Coil inside a permanent magnet. For the most accurate weighings
	<b>Analogue interface:</b> to connect a suitable peripheral device for analogue processing of the measurements		<b>Percentage determination:</b> Determining the deviation in % from the target value (100 %)		<b>Weighing principle: Single cell technology:</b> Advanced version of the force compensation principle with the highest level of precision
	<b>Interface for second balance:</b> For direct connection of a second balance		<b>Weighing units:</b> Can be switched to e.g. nonmetric units at the touch of a key. See balance model. Please refer to KERN's website for more details		<b>Verification possible:</b> The time required for verification is specified in the pictogram
	<b>Network interface:</b> For connecting the scale to an Ethernet network		<b>Weighing with tolerance range:</b> (Checkweighing) Upper and lower limiting can be programmed individually, e.g. for sorting and dosing. The process is supported by an audible or visual signal, see the relevant model		<b>DAKkS calibration possible:</b> The time required for DAKkS calibration is shown in days in the pictogram
	<b>Wireless data transfer:</b> between the weighing unit and the evaluation unit using an integrated radio module		<b>Hold function:</b> (Animal weighing program) When the weighing conditions are unstable, a stable weight is calculated as an average value		<b>Package shipment:</b> The time required for internal shipping preparations is shown in days in the pictogram
					<b>Pallet shipment:</b> The time required for internal shipping preparations is shown in days in the pictogram

\*The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by KERN & SOHN GmbH is under license. Other trademarks and trade names are those of their respective owners.

## KERN – Precision is our business

To ensure the high precision of your balance KERN offers you the the appropriate test weight in the international OIML error limit classes E1-M3 from 1 mg - 2500 kg. In combination with a DAKkS calibration certificate the best pre-requisite for proper balance calibration.

The KERN DAKkS calibration laboratory today is one of the most modern and best-equipped DAKkS calibration laboratories for balances, test weights and force-measurement in Europe.

Thanks to the high level of automation, we can carry out DAKkS calibration of balances, test weights and force-measuring devices 24 hours a day, 7 days a week.

## Range of services:

- DAKkS calibration of balances with a maximum load of up to 50 t
- DAKkS calibration of weights in the range of 1 mg - 2500 kg
- Volume determination and measuring of magnetic susceptibility (magnetic characteristics) for test weights
- Database supported management of checking equipment and reminder service
- Calibration of force-measuring devices
- DAKkS calibration certificates in the following languages DE, GB, FR, IT, ES, NL, PL
- Conformity evaluation and reverification of balances and test weights

## Your KERN specialist dealer:



# WolfLabs

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

**[www.wolflabs.co.uk](http://www.wolflabs.co.uk)**

**Tel : 01759 301142**

**Fax : 01759 301143**

**[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)**

Please contact us if this literature doesn't answer all your questions.