



### **Built for Precision, Powered by Innovation**

The new Explorer (EXR series) offers an intuitive interface, outstanding weighing performance, and an environmentally friendly design. The new generation Explorer establishes a groundbreaking standard for laboratory operations.

### **Unique Features:**

### High Performance Weighing

Achieve consistent results with built-in automatic calibration, touchless sensors and a ventilated weighing pan. Get reliable data with the digital high-speed load cell.

#### Refined User Experience

An intuitive 7-inch glass touchscreen, digital leveling, and simplified data management is easy to use. With an effortless installation and chemical-protective design it ensures efficiency and comfort.

#### Effective Communication Options

Supports HID connection to PC or keyboard without drivers, and provides 2-way communication with RS232, USB, Ethernet, Bluetooth, and WiFi for seamless data transfer.

# **EXPLORER**<sup>™</sup> Semi-Micro Balances


## **Specifications**



Model	EXR125D	EXR125	EXR225D	EXR225
	EXR125DM	EXR125M	EXR225DM	EXR225M
Capacity	82 g / 120 g	120 g	120 g / 220 g	220 g
Readability d, Fine Range	0.01 mg	0.01 mg	0.01 mg	0.01 mg
Readability d, Full Range	0.1 mg	0.01 mg	0.1 mg	0.01 mg
Repeatability (sd.), ≤5% of Full Load Typical	0.01 mg	0.01 mg	0.01 mg	0.01 mg
Repeatability (sd.), 5% of Full Load to Fine Range Max	0.02 mg	0.02 mg	0.02 mg	0.02 mg
Repeatability (sd. ), Fine Range Max to Full Range	0.1 mg	0.02 mg	0.1 mg	0.1 mg
Linearity Deviation Typical	0.06 mg	0.06 mg	0.06 mg	0.06 mg
Linearity Deviation	0.2 mg	0.1 mg	0.2 mg	0.1 mg
Verification Interval, e (EXPM models only)	1 mg	1 mg	1 mg	1 mg
Approval Class (EXPM models only)	Class I	Class I	Class I	Class I
Span Calibration Points	25g, 50g	, 75g, 100g	50g, 100g,	150g, 200g
Weighing Units	g, mg, ct, N, oz, ozt, Grain, dwt, mo, msg, tcl, tola, baht, Custom 1, Custom 2			
Weighing Units (EXRM model only)	g, mg, ct			
Applications	Basic Weighing, Parts Counting, Check Counting, Percent Weighing, Check Weighing, Dynamic Weighing, Totalization, Formulation, Differential, Density Determination, Peak Hold			
Stabilization Time (typical)	0.1 mg: ≤2 seconds, 0.01 mg: ≤ 5 seconds			
Sensitivity Temperature Drift (ppm/°C)	0.8			
Min-Weight (typical) (K=2, U=1%)	2 mg			
Min-Weight (typical) (USP, K=2, U=0.10%)	20 mg			
Min-Weight (optimal) (USP, K=2, U=0.10%, SRP≤0.41d)	8.2 mg			
Terminal Display	7-inch TFT 16.7M Color Touch Screen with Glass Overlay, Detachable Terminal			
Display Resolution	800 x 480 Pixels			
Communication	USB Host (Type A) x 2, USB Device (Type C) x 1 Ethernet (RJ45) x 1, RS232 x 1; Optional Wi-Fi, Bluetooth Dongle			
Leveling System	2 Leveling Feet with Digital Leveling Bubble			
Power Supply	AC Adapter Input: 100-240 VAC 0.5A 50-60 Hz AC Adapter Output: 12 VDC 1.5A			
Weighing Pan Size	Aero Pan 80mm x 80mm			
Touchless Sensor	2 Touchless Sensors on the terminal			
Draftshield with anti-static coating	Standard			
Terminal Status Lights	Standard, 7-color Selectable			
Assembled Dimensions (W x D x H)	230mm x 393mm x 370mm			
Shipping Dimensions (W x D x H)	415mm x 630mm x 561mm			
Net Weight	7.4 kg			
Shipping Weight		10.2	2 kg	

## **EXPLORER**<sup>™</sup> Semi-Micro Balances

### **Outline Dimensions**

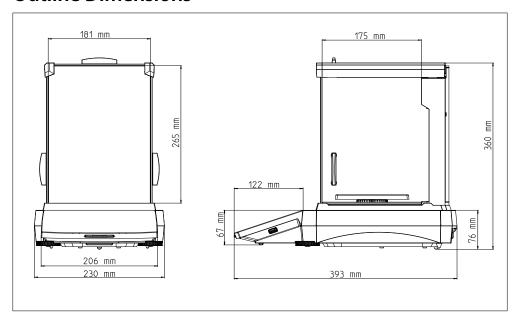

## **Approvals**

- Metrology: OIML R76, EN 45501
- Product Safety: IEC/EN 61010-1; CAN/CSA C22.2 61010-1; UL 61010-1
- Electromagnetic Compatibility: IEC/EN 61326-1 Class B, Basic Environments; FCC Part 15 Class A; Canada ICES-003 Class A

## **EXPLORER**<sup>™</sup> Semi-Micro Balances

### Accessories

a a	Item No.	30944359			
	Accessory Name	Bluetooth Printer, Impact, SF40A/BT			
	Item No.	30064202			
	Accessory Name	Printer, Impact, SF40A			
W or	Item No.	31059239			
	Accessory Name	BT dongle and Wi-Fi dongle LM842			
1000	Item No.	30095929			
	Accessory Name	Static Ionizer, ION-100A			
	Item No.	31059237			
	Accessory Name	Weighing Kit			
	Item No.	80253384			
	Accessory Name	Density Kit, Solids			
	Item No.	83034024			
	Accessory Name	Sinker, Glass, Liquid, Density Kit			
J. H	Item No.	12120799			
	Accessory Name	SF40A Paper roll (57,5mm 2pcs)			
	Item No.	12120798			
	Accessory Name	SF40A Ink Ribbon Cassette			





Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

## www.wolflabs.co.uk

Tel: 01759 301142

Fax: 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.