




MULTI-MODE ROTATOR GENIE™

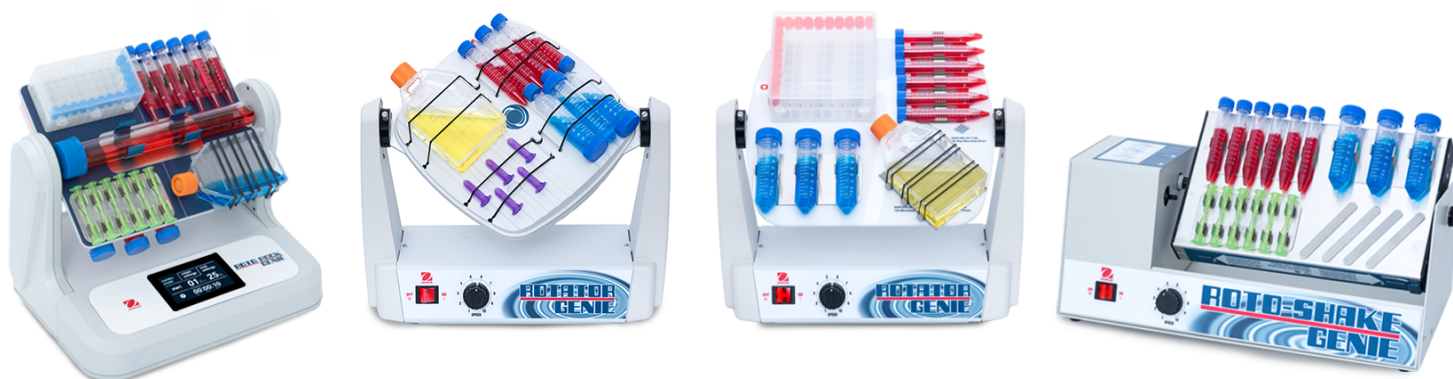

Magneto-magic and Rugged Lab Rotator to Rotate, Rock, or 3D Orbit any Container

Genie rotators combine rotating and rocking in a compact, versatile design. The magneto-magic system features a magnetic platform, clip plates, and bag strips that firmly hold various containers rotating at any angle, in any combination. Rocking tray fits large open containers. 3D and horizontal orbital attachments turn Roto-Shake Genie into an orbital shaker. Touchscreen control allows Roto-Rock Genie to program combined rock and rotate cycles.

Unique Features:

- **Unique Magneto-magic System Offers Unlimited Combinations and Flexibility**

Multi-spec clip plates and strips hold a wide variety of tubes, flasks, microplates, bags, and odd-shaped vessels. Freely attach to the magnetic platform at any angle, in any combination.

- **Mix Smarter in Tight Spaces**

Save space with a compact unit offering multiple mixing modes. Rotator Genie features switchable platforms. Standard platform adjusts to various rotation angles. Usable in cold rooms and incubators.

- **Touchscreen Advantages of Roto-Rock Genie**

Easily switches between rotating and rocking via the touchscreen, and allows programming of speed, time, rocking angle, and combination rock/rotate cycles.

MULTI-MODE ROTATOR GENIE™

Specifications

Model	Roto-Rock Genie	Rotator Genie	Roto-Shake Genie
Control	Touch Control	Analog	
Display	LCD	-	
Accessories Included	Magnetic Platform, Rocking Tray	Standard Platform with Retaining Bands or Magnetic Platform	Clip Plate for 10-13 mm tubes, Clip Plate for 15-17 mm tubes, Clip Plate for 28-30 mm tubes, 4x Bag Mounting Strips, Non-Skid Tray
Capacity	4.5 kg centered on axis for rotating / rocking	2.25 kg centered on axis	4.5 kg centered on axis for rotating / rocking
Motion	Rock, Rotate, Timed auto-switch Rock/Rotate	Rotate	Rock, Rotate, Shake with Accessories
Rotating Speed	1 rpm - 35 rpm	5 rpm - 35 rpm	5 rpm - 35 rpm
Rocking Speed	1 rpm - 70 rpm	-	10 rpm - 70 rpm
Tilt Angle	0 - 25°	-	10°
Timer	1 second - 99 hours, 59 minutes	-	
Magnetic Platform Size (L x W)	250 x 320 mm	290 x 315 mm	200 x 340 mm
Tray Size (L x W)	300 x 300 mm	260 x 260 mm	305 x 305 mm
Safety Certification	CE		
Power	24 V DC, 2.5A	120V, 0.5A, 60Hz 230V, 0.25A, 50Hz 100V, 0.5A, 50/60Hz	120V, 60Hz 230V, 50Hz 100V, 50/60Hz
Dimensions (L x W x H)	300 x 400 x 380 mm	175 x 400 x 400 mm	252 x 508 x 267 mm
Net Weight	9.7 kg	5.2 kg	11.8 kg
Working Environment	0°C - 38°C, 95%RH, non-condensing		

Accessories



83041252 (REF #: SI-R150)

Stacking Tray Assembly, Roto-Rock

Compatibility: Roto-Rock Genie



83041367 (REF #: SI-1175)

Horizontal Orbital Shaker Attachment

Compatibility: Roto-Shake Genie



83041368 (REF #: SI-1170)

Stackable Tray for Roto-Shake Genie

Compatibility: Roto-Shake Genie



83041369 (REF #: SI-1150)

3-D Orbital Shaker Attachment


Compatibility: Roto-Shake Genie



83041370 (REF #: SI-1140)

Magnetic 50 Microtube Box


Compatibility: Roto-Rock Genie, Rotator Genie*

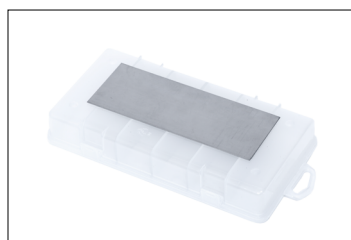

83041371 (REF #: SI-1139)

Magnetic Strip Blot Box

Compatibility: Roto-Rock Genie

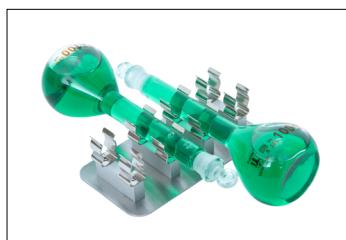

83041372 (REF #: SI-1136)

1 Clip Plate for 100ml Volumetric Flasks

Compatibility: Roto-Rock Genie, Rotator Genie*, Roto-Shake Genie




83041373 (REF #: SI-1135)

1 Magnetic Covered 100 Microtube Box

Compatibility: Roto-Rock Genie, Rotator Genie*, Roto-Shake Genie


MULTI-MODE ROTATOR GENIE™

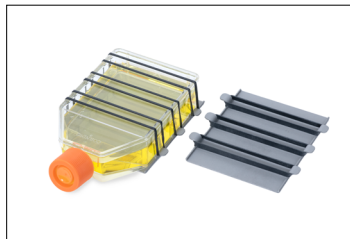

83041374 (REF #: SI-1134)

2 Universal Clip Plates

Compatibility: Roto-Rock Genie, Rotator Genie*, Roto-Shake Genie

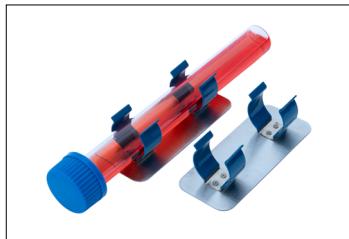

83041377 (REF #: SI-1130)

2 Clip Plates for Hybridization Tubes


Compatibility: Roto-Rock Genie, Rotator Genie*, Roto-Shake Genie



83041378 (REF #: SI-1122)

6 Clip Plates for 3 each 28-30mm tubes


Compatibility: Roto-Rock Genie, Rotator Genie*, Roto-Shake Genie

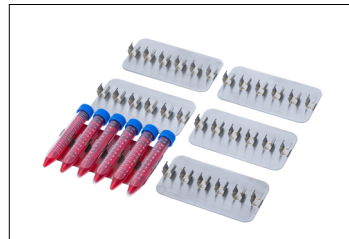

83041379 (REF #: SI-1121)

6 Clip Plates for 6 each 15-17mm tubes


Compatibility: Roto-Rock Genie, Rotator Genie*, Roto-Shake Genie

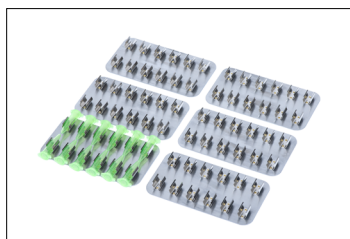

83041380 (REF #: SI-1120)

6 Clip Plates for 12 each 10-13mm tubes

Compatibility: Roto-Rock Genie, Rotator Genie*, Roto-Shake Genie

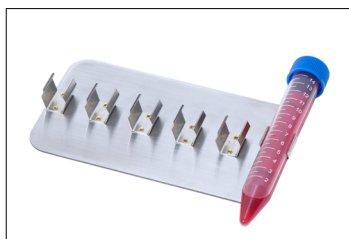

83041443 (REF #: SI-1124)

Clip Plate (1) For 6 Each 15-17mm Tubes

Compatibility: Roto-Rock Genie, Rotator Genie*, Roto-Shake Genie



83041445 (REF #: SI-1125)

Clip Plate (1) For 3 Each 28-30mm Tubes


Compatibility: Roto-Rock Genie, Rotator Genie*, Roto-Shake Genie



83041457 (REF #: SI-2210)

Magnetic Platform for Rotator-Genie


Compatibility: Rotator Genie

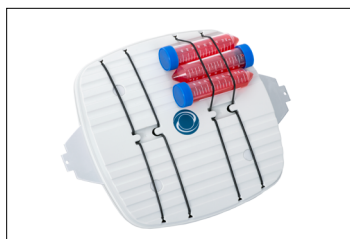

83041458 (REF #: SI-2110)

Standard Platform for Rotator-Genie


Compatibility: Rotator Genie

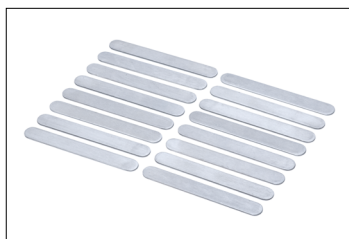

83041465 (REF #: SI-1126)

Bag Mounting Strips (pack of 16)

Compatibility: Roto-Rock Genie, Rotator Genie*, Roto-Shake Genie

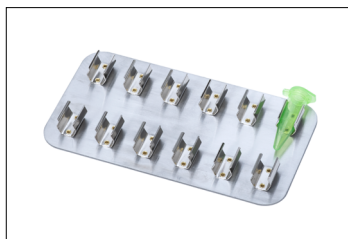

83041466 (REF #: SI-1123)

1 Clip Plate for 12 each 10-13mm tubes

Compatibility: Roto-Rock Genie, Rotator Genie*, Roto-Shake Genie



83041723 (REF #: 0A-1100-070)

Non-Skid Tray

Compatibility: Roto-Shake Genie

* Compatible with Rotator Genie Magnetic Platform, not compatible with Rotator Genie Standard Platform

Roto-Rock Genie Use with Accessory Attachments



with included Rocking Tray



with 83041252 (SI-R150)



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.